

# INTERNATIONAL STANDARD

**IEC**  
**61223-2-11**

First edition  
1999-09

---

---

## **Evaluation and routine testing in medical imaging departments –**

### **Part 2-11: Constancy tests – Equipment for general direct radiography**

*Essais d'évaluation et de routine dans  
les services d'imagerie médicale –*

*Partie 2-11:  
Essais de constance –  
Appareils de radiographie générale directe*



Reference number  
IEC 61223-2-11:1999(E)

## Numbering

As from 1 January 1997 all IEC publications are issued with a designation in the 60000 series.

## Consolidated publications

Consolidated versions of some IEC publications including amendments are available. For example, edition numbers 1.0, 1.1 and 1.2 refer, respectively, to the base publication, the base publication incorporating amendment 1 and the base publication incorporating amendments 1 and 2.

## Validity of this publication

The technical content of IEC publications is kept under constant review by the IEC, thus ensuring that the content reflects current technology.

Information relating to the date of the reconfirmation of the publication is available in the IEC catalogue.

Information on the subjects under consideration and work in progress undertaken by the technical committee which has prepared this publication, as well as the list of publications issued, is to be found at the following IEC sources:

- **IEC web site\***
- **Catalogue of IEC publications**  
Published yearly with regular updates  
(On-line catalogue)\*
- **IEC Bulletin**  
Available both at the IEC web site\* and as a printed periodical

## Terminology, graphical and letter symbols

For general terminology, readers are referred to IEC 60050: *International Electrotechnical Vocabulary* (IEV).

For graphical symbols, and letter symbols and signs approved by the IEC for general use, readers are referred to publications IEC 60027: *Letter symbols to be used in electrical technology*, IEC 60417: *Graphical symbols for use on equipment*, *Index, survey and compilation of the single sheets* and IEC 60617: *Graphical symbols for diagrams*.

\* See web site address on title page.

# INTERNATIONAL STANDARD

**IEC**  
**61223-2-11**

First edition  
1999-09

---

---

## **Evaluation and routine testing in medical imaging departments –**

### **Part 2-11: Constancy tests – Equipment for general direct radiography**

*Essais d'évaluation et de routine dans  
les services d'imagerie médicale –*

*Partie 2-11:  
Essais de constance –  
Appareils de radiographie générale directe*

© IEC 1999 — Copyright - all rights reserved

No part of this publication may be reproduced or utilized in any form or by any means, electronic or mechanical, including photocopying and microfilm, without permission in writing from the publisher.

International Electrotechnical Commission  
Telefax: +41 22 919 0300

3, rue de Varembé Geneva, Switzerland  
e-mail: [inmail@iec.ch](mailto:inmail@iec.ch) IEC web site <http://www.iec.ch>



Commission Electrotechnique Internationale  
International Electrotechnical Commission  
Международная Электротехническая Комиссия

PRICE CODE

**W**

*For price, see current catalogue*

# CONTENTS

|                                                                                | Page |
|--------------------------------------------------------------------------------|------|
| FOREWORD .....                                                                 | 3    |
| Clause                                                                         |      |
| 1 Scope and object .....                                                       | 5    |
| 1.1 Scope .....                                                                | 5    |
| 1.2 Object .....                                                               | 5    |
| 2 Normative references .....                                                   | 6    |
| 3 Terminology .....                                                            | 6    |
| 3.1 Degree of requirements .....                                               | 6    |
| 3.2 Use of terms .....                                                         | 7    |
| 4 General aspects of CONSTANCY TESTS .....                                     | 7    |
| 4.1 General conditions affecting test procedures .....                         | 8    |
| 4.2 Establishment of BASELINE VALUES .....                                     | 8    |
| 4.3 Frequency of CONSTANCY TESTS .....                                         | 8    |
| 4.4 Identification of equipment, instrumentation and test conditions .....     | 8    |
| 4.5 Measured functional parameters .....                                       | 9    |
| 5 Performance tests .....                                                      | 10   |
| 5.1 RADIATION output from the X-RAY SOURCE ASSEMBLY .....                      | 10   |
| 5.2 RADIATION input to the IMAGE RECEPTION AREA .....                          | 12   |
| 5.3 Geometric characteristics .....                                            | 14   |
| 5.4 Resolution of high-contrast detail .....                                   | 19   |
| 5.5 Variation in optical density throughout a RADIOGRAM .....                  | 20   |
| 6 Statement of compliance .....                                                | 22   |
| Figure 1 Example of a film marker TEST DEVICE .....                            | 23   |
| Figure 2 Example of a TEST DEVICE for perpendicular position .....             | 24   |
| Figure 3 Example of an alignment TEST DEVICE .....                             | 25   |
| Figure 4 Examples of a high-contrast TEST DEVICE .....                         | 26   |
| Figure 5 Arrangements for testing geometric conditions .....                   | 28   |
| Figure 6 Geometric coincidences .....                                          | 29   |
| Annex A (normative) Terminology – Index of defined terms .....                 | 30   |
| Annex B (informative) Example of a form for the standardized test report ..... | 32   |
| Annex C (informative) Guidance on action to be taken .....                     | 35   |
| Annex D (normative) PHANTOMS and TEST DEVICES .....                            | 36   |

## INTERNATIONAL ELECTROTECHNICAL COMMISSION

## EVALUATION AND ROUTINE TESTING IN MEDICAL IMAGING DEPARTMENTS –

**Part 2-11: Constancy tests –  
Equipment for general direct radiography**

## FOREWORD

- 1) The IEC (International Electrotechnical Commission) is a worldwide organization for standardization comprising all national electrotechnical committees (IEC National Committees). The object of the IEC is to promote international co-operation on all questions concerning standardization in the electrical and electronic fields. To this end and in addition to other activities, the IEC publishes International Standards. Their preparation is entrusted to technical committees; any IEC National Committee interested in the subject dealt with may participate in this preparatory work. International, governmental and non-governmental organizations liaising with the IEC also participate in this preparation. The IEC collaborates closely with the International Organization for Standardization (ISO) in accordance with conditions determined by agreement between the two organizations.
- 2) The formal decisions or agreements of the IEC on technical matters express, as nearly as possible, an international consensus of opinion on the relevant subjects since each technical committee has representation from all interested National Committees.
- 3) The documents produced have the form of recommendations for international use and are published in the form of standards, technical reports or guides and they are accepted by the National Committees in that sense.
- 4) In order to promote international unification, IEC National Committees undertake to apply IEC International Standards transparently to the maximum extent possible in their national and regional standards. Any divergence between the IEC Standard and the corresponding national or regional standard shall be clearly indicated in the latter.
- 5) The IEC provides no marking procedure to indicate its approval and cannot be rendered responsible for any equipment declared to be in conformity with one of its standards.
- 6) Attention is drawn to the possibility that some of the elements of this International Standard may be the subject of patent rights. The IEC shall not be held responsible for identifying any or all such patent rights.

International Standard IEC 61223-2-11 has been prepared by subcommittee 62B: Diagnostic imaging equipment, of IEC technical committee 62: Electrical equipment in medical practice.

The text of this standard is based on the following documents:

| FDIS         | Report on voting |
|--------------|------------------|
| 62B/373/FDIS | 62B/385/RVD      |

Full information on the voting for the approval of this standard can be found in the report on voting indicated in the above table.

This publication has been drafted in accordance with the ISO/IEC Directives, Part 3.

Annexes A and D form an integral part of this standard.

Annexes B and C are for information only.

This standard forms part 2-11 of IEC 61223, which will include the following parts:

- Part 1: General aspects
- Part 2-1: Constancy tests – Film processors
- Part 2-2: Constancy tests – Radiographic cassettes and film changers – Film-screen contact and relative sensitivity of the screen-cassette assembly
- Part 2-3: Constancy tests – Darkroom safelight conditions
- Part 2-4: Constancy tests – Hard copy cameras
- Part 2-5: Constancy tests – Image display devices
- Part 2-6: Constancy tests – X-ray equipment for computed tomography
- Part 2-7: Constancy tests – Equipment for intra-oral dental radiography excluding dental panoramic equipment
- Part 2-9: Constancy tests – Equipment for indirect radioscopy and indirect radiography
- Part 2-10: Constancy tests – X-ray equipment for mammography
- Part 2-11: Constancy tests – Equipment for general direct radiography

The committee has decided that this publication remains valid until 2003. At this date, in accordance with the committee's decision, the publication will be

- reconfirmed;
- withdrawn;
- replaced by a revised edition, or
- amended.

A bilingual version of this standard may be issued at a later date.

**EVALUATION AND ROUTINE TESTING IN MEDICAL IMAGING DEPARTMENTS –****Part 2-11: Constancy tests –  
Equipment for general direct radiography****1 Scope and object****1.1 Scope**

This part of IEC 61223 applies to those components of X-RAY EQUIPMENT which

- generate, influence the propagation of, and detect X-RADIATION; and
- process, present and store radiographic information in RADIOLOGICAL INSTALLATIONS with diagnostic X-ray systems using RADIOGRAPHIC FILM in DIRECT RADIOGRAPHY.

This standard is a part of a series of Particular Publications (International standards and technical reports), which define methods of testing the constancy of operation of various subsystems of diagnostic X-RAY EQUIPMENT.

This standard does not apply to equipment for special applications such as mammographic X-RAY EQUIPMENT or dental X-RAY EQUIPMENT; see complete list of all parts 2 of IEC 61223 in the foreword.

This standard gives methods of tests for the constancy of properties of diagnostic X-RAY EQUIPMENT as described in IEC 61223-1 (see clause 2).

This part of IEC 61223 is designed to be applicable to equipment for general direct radiography without digital imaging devices.

**1.2 Object**

This standard defines

- the essential parameters which describe or affect the performance of the above components of X-RAY EQUIPMENT;
- methods of checking that variations in measured quantities related to those parameters are within acceptable limits, in order to maintain adequate standards of imaging whilst reducing unnecessary IRRADIATION of the PATIENT.

The methods are based upon assessments of RADIOGRAMS of appropriate TEST DEVICES.

The purpose of the methods is

- to establish a reference level of performance when such equipment is accepted;
- to detect and verify any significant variation in performance which may require corrective action.

Because RADIOLOGICAL INSTALLATIONS differ widely from each other, it is not possible in this standard to specify target values and tolerances for the parameters which would be generally applicable as criteria of acceptable performance. Guidance is given, however, as to the degree of variation in single measurements which might require appropriate action.

This standard does not deal with

- aspects of mechanical and electrical safety;
- checks of the effectiveness of the direct means of protection against X-RADIATION;
- optimization of imaging performance.

With regard to the measurements, reference is made to methods described in related publications, which for practical reasons should be carried out prior to the application of the methods described in this standard (see clause 2).

## 2 Normative references

The following normative documents contain provisions which, through reference in this text, constitute provisions of this International Standard. For dated references, subsequent amendments to, or revisions of, any of these publications do not apply. However, parties to agreements based on this International Standard are encouraged to investigate the possibility of applying the most recent editions of the normative documents indicated below. For undated references, the latest edition of the normative document referred to applies. Members of IEC and ISO maintain registers of currently valid International Standards.

IEC 60601-1-3:1994, *Medical electrical equipment – Part 1: General requirements for safety*  
3. Collateral standard: *General requirements for radiation protection in diagnostic X-ray equipment*

IEC 60788:1984, *Medical radiology – Terminology*

IEC 61223-1:1993, *Evaluation and routine testing in medical imaging departments – Part 1: General aspects*

IEC 61223-2-1:1993, *Evaluation and routine testing in medical imaging departments – Part 2-1: Constancy tests – Film processors*

IEC 61223-2-2:1993, *Evaluation and routine testing in medical imaging departments – Part 2-2: Constancy tests – Radiographic cassettes and film changers – Film-screen contact and relative sensitivity of the screen-cassette assembly*

IEC 61223-2-3:1993, *Evaluation and routine testing in medical imaging departments – Part 2-3: Constancy tests – Darkroom safelight conditions*

## 3 Terminology

### 3.1 Degree of requirements

In this standard, certain terms (which are not printed in SMALL CAPITALS) have particular meanings, as follows:

- "shall" indicates a requirement that is mandatory for compliance;
- "should" indicates a strong recommendation that is not mandatory for compliance;
- "may" indicates a permitted manner of complying with a requirement or of avoiding the need to comply;
- "specific" is used to indicate definitive information stated in this standard or referenced in other standards, usually concerning particular operating conditions, test arrangements or values connected with compliance;

- "specified" is used to indicate definitive information stated by the manufacturer in accompanying documents or in other documentation relating to the equipment under consideration, usually concerning its intended purpose, or the parameters or conditions associated with its use or with testing to determine compliance.

### 3.2 Use of terms

In this standard, terms printed in SMALL CAPITALS are used as defined in IEC 60788 or other IEC publications; see annex A. Where a defined term is used as a qualifier in another defined or undefined term, it is not printed in SMALL CAPITALS, unless the concept thus qualified is defined, or recognized as a "derived term without definition". Test specifications are in italics.

NOTE – Attention is drawn to the fact that, in cases where the concept addressed is not strongly confined to the definition given in one of the publications listed above, a corresponding term is printed in lower-case letters.

## 4 General aspects of CONSTANCY TESTS

The methods for testing the constancy described in this standard are intended to enable the OPERATOR to detect changes in image quality of images produced by X-RAY EQUIPMENT.

For the results of the CONSTANCY TESTS described in this standard to be valid, it is essential to ensure that they are not significantly influenced by anything other than changes in the parameters under test.

In particular, attention shall be paid to darkroom safelight conditions, according to IEC 61223-2-3, and proper film processing, according to IEC 61223-2-1 (see clause 2). When using FILM ILLUMINATORS special attention should also be paid to lighting conditions.

Careful consideration shall be given to the operating and test conditions under which the equipment is checked, including environmental influences.

All equipment under test and the test equipment shall be identified at the initial CONSTANCY TEST in order to ensure that the same items are used in subsequent CONSTANCY TESTS.

Appropriate individual TEST DEVICES are described in detail in annex D for each of the test measurements described in clause 5 of this standard. In practice, a composite TEST DEVICE may be used, which combines the properties of some of the individual TEST DEVICES. Furthermore, some of the test procedures described separately in clause 5 may be performed simultaneously.

NOTE – If the MANUFACTURER provides proposals for the method and frequency of CONSTANCY TESTS in the ACCOMPANYING DOCUMENTS, they should preferably be followed.

#### 4.1 General conditions affecting test procedures

The CONSTANCY TESTS described in this standard have been designed to be easily reproducible, i.e., their results should be affected only by changes in the parameters under investigation. The number of test tools and test equipment has been kept to a minimum and restricted where possible to devices that are passive, inherently simple or reasonably stable. It is important

- to perform CONSTANCY TESTS with LOADING FACTORS which are the same as those used most frequently in clinical practice;
- to record and reproduce all significant settings of the X-RAY EQUIPMENT and ACCESSORIES each time a test is undertaken and to check that the same equipment, components and ACCESSORIES are being used;
- to consider particularly variations in supply voltage, on the results;
- to use RADIOGRAPHIC FILM which is handled, processed and viewed in accordance with the standards and technical reports referenced in clause 2;
- to check the performance of the test instrumentation regularly, particularly when any significant variation in the X-RAY EQUIPMENT is suspected.

NOTE – When appropriate national standards exist, measuring equipment should be referable to them.

Before the CONSTANCY TESTS are started, the constancy of the radiographic cassettes, the radiographic film, the film processing and the film viewing conditions shall be checked.

#### 4.2 Establishment of BASELINE VALUES

When new X-RAY EQUIPMENT is brought into use, or any component of the X-RAY EQUIPMENT, ACCESSORIES or test equipment is changed, which may cause a variation in the test result, an initial CONSTANCY TEST shall be carried out immediately after an ACCEPTANCE TEST or STATUS TEST has indicated that the performance is satisfactory. The purpose of the initial CONSTANCY TEST is to establish new BASELINE VALUES for the parameters tested.

#### 4.3 Frequency of CONSTANCY TESTS

The CONSTANCY TESTS shall be repeated as directed in the appropriate subclauses of this standard. In addition, the CONSTANCY TESTS should be repeated

- whenever malfunction is suspected;
- immediately after the equipment has undergone maintenance that could affect the performance parameter under test;
- to confirm the test results, whenever the results are outside the criteria.

Records of the BASELINE VALUES shall be kept until a new initial CONSTANCY TEST is performed. The results of the CONSTANCY TESTS shall be kept at least two years.

#### 4.4 Identification of equipment, instrumentation and test conditions

All X-RAY EQUIPMENT under test or used for testing shall be unequivocally identified.

Interchangeable components of X-RAY EQUIPMENT such as

- ADDED FILTERS;
- BEAM LIMITING DEVICE;
- PATIENT SUPPORT or other attenuating material in the RADIATION BEAM;
- RADIATION DETECTOR of the AUTOMATIC CONTROL SYSTEM;
- ANTI-SCATTER GRID;

together with items of test instrumentation such as

- combination of RADIOGRAPHIC CASSETTE and INTENSIFYING SCREENS;
- TEST DEVICES;
- RADIOGRAPHIC FILM type and emulsion number;
- FILM PROCESSOR;
- sensitometer;
- densitometer;

and settings of variables such as

- FOCAL SPOT TO IMAGE RECEPTOR DISTANCE;
- AUTOMATIC CONTROL SYSTEM density control and sensor position;
- LOADING FACTORS;
- nominal FOCAL SPOT size, if applicable;

shall be marked and/or recorded so that the items and settings used in the initial CONSTANCY TEST can be used with the equipment under test.

NOTE 1 – Most of the tests should be performed with the particular RADIOGRAPHIC CASSETTE used for the initial CONSTANCY TEST. This cassette, hereafter referred to as the "test cassette", may be kept exclusively for test procedures or identified among those used regularly in clinical work. Whereas the first is more likely to provide a stable tool for revealing changes in the equipment, the second approach will be more representative of changes of the whole system, including those due to the ageing of the cassette itself.

NOTE 2 – It is essential that any RADIOGRAPHIC FILM used in the test is of the same type as the film used for the CONSTANCY TEST for the FILM PROCESSOR.

RADIOGRAMS of a TEST DEVICE shall be taken with the test cassette, with the same INTENSIFYING SCREENS and film type. The RADIOGRAPHIC FILM shall be processed under known comparable conditions and appropriate allowance shall be made for any changes in film batch according to the method specified in IEC 61223-2-1. If the film type and/or associated processing conditions are changed, a new initial CONSTANCY TEST shall be carried out.

#### 4.5 Measured functional parameters

The imaging performance of a direct radiographic system is considered to be constant if the variations of the following functional parameters are found to meet applicable criteria:

- RADIATION output from the X-RAY SOURCE ASSEMBLY; see 5.1;
- RADIATION input to the IMAGE RECEPTION AREA; see 5.2;
- geometric characteristics; see 5.3;
- resolution of high-contrast detail; see 5.4;
- variation in optical density throughout a RADIOGRAM; see 5.5.

The constancy of the imaging performance may be affected by variations in one or more of the following parameters:

- value and wave-form of the mains voltage;
- value and wave-form of the X-RAY TUBE VOLTAGE;
- X-RAY TUBE CURRENT;
- LOADING TIME;
- IRRADIATION TIME under AUTOMATIC EXPOSURE CONTROL;
- FILTRATION and ATTENUATION by layers in the RADIATION BEAM;
- roughening of the TARGET in the ANODE of the X-RAY TUBE;

- distances of planes of interest from the FOCAL SPOT;
- limitation of the RADIATION BEAM;
- direction of the RADIATION BEAM;
- alignment of the RADIATION BEAM with the IMAGE RECEPTION AREA;
- coincidence of the RADIATION FIELD with the field indicated by the LIGHT FIELD-INDICATOR;
- movement of the MOVING GRID;
- positioning of the ANTI-SCATTER GRID;
- imaging properties of the FOCAL SPOT;
- mechanical stability.

## 5 Performance tests

### 5.1 RADIATION output from the X-RAY SOURCE ASSEMBLY

#### 5.1.1 Summary

The RADIATION output from the X-RAY SOURCE ASSEMBLY is measured by means of a RADIATION METER. The measurements are carried out under manual and/or AUTOMATIC EXPOSURE CONTROL, depending upon the type and use of the X-RAY EQUIPMENT.

#### 5.1.2 Test equipment

Measurements are performed using an integrating RADIATION METER having an overall reproducibility of  $\pm 5$  % (including long-term stability, instrument noise and read-out ability).

As a substitute for the PATIENT when testing under AUTOMATIC EXPOSURE CONTROL, the ATTENUATION PHANTOM shall be used. This is intended to provide appropriate ATTENUATION and hardening of the X-RAY BEAM. A detailed description including consideration for the correct choice of an ATTENUATION PHANTOM is provided in annex D.

#### 5.1.3 Test procedure

*Place the RADIATION DETECTOR of the RADIATION METER in the RADIATION BEAM emerging from the X-RAY SOURCE ASSEMBLY. Measurement geometry shall be such that*

- *the distance of the RADIATION DETECTOR from the FOCAL SPOT;*
- *the position of the RADIATION DETECTOR within the RADIATION FIELD;*

*are reproduced to within  $\pm 1$  % of the FOCAL SPOT to RADIATION DETECTOR distance used in the initial CONSTANCY TEST. The same RADIATION FIELD size shall be used.*

*Whenever possible, carry out this test under both manual and AUTOMATIC EXPOSURE CONTROL. If a measurement according to 5.2.3.2 is performed, a measurement according to 5.1.3.2 may be regarded as redundant.*

##### 5.1.3.1 Testing under manual control

*Operate the X-RAY EQUIPMENT using manual settings of the LOADING FACTORS identical to those used in the initial CONSTANCY TEST.*

*Record the reading of the RADIATION METER.*

### 5.1.3.2 Testing under AUTOMATIC EXPOSURE CONTROL

*Align the X-RAY SOURCE ASSEMBLY with the X-RAY IMAGE RECEPTOR as in normal clinical practice.*

*Place the ATTENUATION PHANTOM in the RADIATION BEAM between the RADIATION DETECTOR of the RADIATION METER and the RADIATION DETECTORS of the AUTOMATIC CONTROL SYSTEM.*

*Position the RADIATION DETECTOR of the RADIATION METER so as not to affect the operation of the AUTOMATIC CONTROL SYSTEM.*

*NOTE – For every RADIATION DETECTOR, a specific instruction has to be given as to where it has to be positioned in order to minimize the effect on the AUTOMATIC EXPOSURE CONTROL. If ATTENUATION PHANTOMS are used which are attached to the X-RAY SOURCE ASSEMBLY (A2 or A3), making the measurement behind the PHANTOM (if the sensitivity of the RADIATION DETECTOR allows it) has to be considered.*

*Place an unloaded RADIOGRAPHIC CASSETTE in the CASSETTE CHANGER and operate the X-RAY EQUIPMENT in conjunction with the AUTOMATIC CONTROL SYSTEM using settings of the X-RAY TUBE VOLTAGE identical to those used in the initial CONSTANCY TEST. Use the same RADIOGRAPHIC CASSETTE.*

*Record the reading of the RADIATION METER. In addition, record the IRRADIATION TIME, CURRENT TIME PRODUCT, etc. after each IRRADIATION has been made, if these are indicated.*

### 5.1.4 Data evaluation

*Compare the measured values of RADIATION output with the established BASELINE VALUES.*

### 5.1.5 Criteria to be applied

#### 5.1.5.1 Testing under manual control

The RADIATION output should be within  $\pm 20$  % of the BASELINE VALUES.

#### 5.1.5.2 Testing under AUTOMATIC EXPOSURE CONTROL

The criteria to be applied depend on material used in the ATTENUATION PHANTOM.

If materials with low atomic number (up to 14) are used (for example water, polymethyl-methacrylate (PMMA), aluminium) the RADIATION output should be within the range +25 % to –20 % of the BASELINE VALUES.

For high atomic number materials, for example copper or lead, the action levels for the RADIATION output should be within  $\pm 25$  % of the BASELINE VALUES.

When lead is used with X-RAY TUBE VOLTAGES above 90 kV, the action levels appropriate to low atomic number materials shall be applied.

### 5.1.6 Action to be taken

If the system fails to meet the criteria, the guidance given in annex C should be followed.

*NOTE – A gradual decrease in the RADIATION output is to be expected as a result of ageing of the X-RAY TUBE. In order to allow for this, it will be necessary from time to time to determine new BASELINE VALUES. This holds especially for testing under manual control, whereas for AUTOMATIC EXPOSURE CONTROL a decrease in the RADIATION output is compensated and cannot be detected.*

### 5.1.7 Frequency of CONSTANCY TESTS

Initially, a series of CONSTANCY TESTS should be carried out daily for at least one week in order to establish the BASELINE VALUES by calculating the mean value of the output measurements.

Subsequently, the CONSTANCY TESTS should be repeated in a two-week cycle up to six months, in order to obtain data regarding the reliability of the X-RAY SOURCE ASSEMBLY, HIGH-VOLTAGE GENERATOR and AUTOMATIC CONTROL SYSTEM. Thereafter, the CONSTANCY TESTS shall be repeated according to the INSTRUCTIONS FOR USE as provided by the MANUFACTURER. If no such information is provided, the CONSTANCY TESTS shall be performed at least quarterly.

## 5.2 RADIATION input to the IMAGE RECEPTION AREA

### 5.2.1 Summary

The RADIATION input to the IMAGE RECEPTION AREA is determined from measurements of the optical density, recorded at specific points of a RADIOGRAM. Alternatively, the RADIATION input may be directly measured by means of a RADIATION METER.

NOTE – However a RADIATION METER takes no account of the X-RAY BEAM quality. The film method should therefore be preferred.

The measurements are carried out under manual and/or AUTOMATIC EXPOSURE CONTROL, depending upon the type and use of the X-RAY EQUIPMENT.

### 5.2.2 Test equipment

This test shall be performed using the same ANTI-SCATTER GRID, the same dedicated RADIOGRAPHIC CASSETTE and associated screen-film combination of the type normally used with the X-RAY EQUIPMENT under test. The systematic error of the optical densitometer used shall not be higher than  $\pm 0,02$ .

If a RADIATION METER is used, the instrument shall be of the integrating type with an overall reproducibility of  $\pm 5\%$  (including long-term stability, instrument noise and read-out ability).

The ATTENUATION PHANTOM shall be used as a substitute for the PATIENT. This is intended to provide appropriate ATTENUATION and hardening of the X-RAY BEAM. In addition, the film marker TEST DEVICE shall also be used in order to enable specific points in the RADIOGRAM to be identified where measurements of the optical density are to be made.

Detailed descriptions including consideration for the correct choice of the ATTENUATION PHANTOM and the film marker TEST DEVICE are provided in annex D.

### 5.2.3 Test procedure

*If the equipment is provided with AUTOMATIC EXPOSURE CONTROL, perform the CONSTANCY TEST in automatic mode. Whenever possible, repeat this CONSTANCY TEST under manual control.*

*Align the X-RAY SOURCE ASSEMBLY with the X-RAY IMAGE RECEPTOR as in normal clinical practice.*

*If a RADIOGRAPHIC CASSETTE is used, put it into its place. If a RADIATION METER is used, place its RADIATION DETECTOR in the IMAGE RECEPTION AREA. Place the ATTENUATION PHANTOM and the film marker TEST DEVICE in the X-RAY BEAM between the FOCAL SPOT and the X-RAY IMAGE RECEPTOR.*

*For any item of test equipment, the measurement geometry shall be such that the distance from the FOCAL SPOT and the position within the RADIATION FIELD are reproduced to within 1 % of the distance from the FOCAL SPOT used in the initial CONSTANCY TEST. The same size of RADIATION FIELD shall be used.*

*Select at least two X-RAY TUBE VOLTAGE settings that are identical to the settings used in the initial CONSTANCY TEST.*

*a) Testing under manual control*

*Operate the X-RAY EQUIPMENT by using the manual settings of the LOADING FACTORS identical to those in the initial CONSTANCY TEST.*

*Process the irradiated film in accordance with the procedure referred to in 4.4.*

*Measure the optical density at specific points of the RADIOGRAM or record the signal indicated by the RADIATION METER.*

*If necessary, adjust the measured values of density in accordance with the procedure referred to in 4.4 to allow for any changes in film batch or processing conditions.*

*b) Testing under AUTOMATIC EXPOSURE CONTROL*

*Operate the X-RAY EQUIPMENT in conjunction with the AUTOMATIC CONTROL SYSTEM and using the X-RAY TUBE VOLTAGE and all other appropriate settings identical to those used in the initial CONSTANCY TEST.*

*Follow the same procedure as indicated in a) above for testing under manual control.*

#### **5.2.4 Data evaluation**

*Compare measured values of optical density or RADIATION input with the established BASELINE VALUES.*

#### **5.2.5 Criteria to be applied**

*a) Testing under manual control*

The optical density should be within  $\pm 0,3$  of the BASELINE VALUE. If a RADIATION METER is used, the RADIATION input should be within  $\pm 30$  % of the BASELINE VALUE.

*b) Testing under AUTOMATIC EXPOSURE CONTROL*

The optical density should be within  $\pm 0,15$  of the BASELINE VALUE. If a RADIATION METER is used, the RADIATION input should be within  $\pm 15$  % of the BASELINE VALUE.

The above tolerances apply to RADIOGRAPHIC FILMS with an average gradient between 2 and 3.

#### **5.2.6 Action to be taken**

If the system fails to meet the criteria, the guidance given in annex C should be followed.

#### **5.2.7 Frequency of CONSTANCY TESTS**

Initially, a series of CONSTANCY TESTS should be carried out daily for at least one week in order to establish the BASELINE VALUES by calculating the mean of the measurements of optical densities or RADIATION input.

Thereafter, the CONSTANCY TESTS shall be repeated according to the INSTRUCTIONS FOR USE as provided by the MANUFACTURER. If no such information is provided, the CONSTANCY TESTS shall be performed at least quarterly.

### 5.3 Geometric characteristics

#### 5.3.1 Summary

These tests shall be performed to check the constancy of the following geometric characteristics of the X-RAY EQUIPMENT:

- the indicated FOCAL SPOT TO IMAGE RECEPTOR DISTANCE;
- the perpendicular position of the RADIATION BEAM AXIS to the IMAGE RECEPTOR PLANE;

NOTE – Generally, it is sufficient to test the REFERENCE DIRECTION of the X-RAY SOURCE ASSEMBLY with respect to the surface of the PATIENT SUPPORT or the ENTRANCE SURFACE of the X-RAY IMAGE RECEPTOR. A change in parallelism between these surfaces and the IMAGE RECEPTOR PLANE is in most cases unlikely.

- the coincidence of the RADIATION FIELD with the LIGHT FIELD; see figure 202 of IEC 60601-1-3;
- the coincidence of the RADIATION FIELD with the X-RAY IMAGE RECEPTOR; see figure 202 of IEC 60601-1-3;
- the read-out of the RADIATION FIELD size, where applicable.

#### 5.3.2 Test equipment

The following test equipment is required:

- tape measure;
- two cassettes with screens, of different sizes, for example 24 cm × 30 cm and 35 cm × 43 cm;
- RADIOGRAPHIC FILM;
- ruler;
- spirit level;
- TEST DEVICE for perpendicular position;
- TEST DEVICE for alignment.

In order to test the perpendicular position of the RADIATION BEAM AXIS to the IMAGE RECEPTOR PLANE, the TEST DEVICE shall be used.

The alignment TEST DEVICE is used to enable the edges and centre of the LIGHT FIELD and the unambiguous orientation of this TEST DEVICE to be identified on the RADIOGRAM.

If IRRADIATIONS are made under AUTOMATIC EXPOSURE CONTROL, the ATTENUATION PHANTOM shall also be used as a substitute for the PATIENT in order to provide appropriate ATTENUATION and hardening of the X-RAY BEAM.

Detailed descriptions of the ATTENUATION PHANTOM and the TEST DEVICES for the perpendicular position and alignment are provided in annex D.

### 5.3.3 Test procedure

#### 5.3.3.1 Indicated FOCAL SPOT TO IMAGE RECEPTOR DISTANCE

- Step 1: Using the indicator associated with the X-RAY EQUIPMENT, set the FOCAL SPOT TO IMAGE RECEPTOR DISTANCE to the value used in the initial CONSTANCY TEST.*
- Step 2: Measure the FOCAL SPOT TO IMAGE RECEPTOR DISTANCE using the tape measure. If the position of the FOCAL SPOT is not clearly indicated on the surface of the X-RAY SOURCE ASSEMBLY, make this measurement from the same point on the X-RAY SOURCE ASSEMBLY as that used in the initial CONSTANCY TEST, for example the nearest point on the X-RAY SOURCE ASSEMBLY to the X-RAY IMAGE RECEPTOR.*

#### 5.3.3.2 Other geometric characteristics of over-couch X-RAY TUBE ASSEMBLIES with manually adjusted, fixed or automatic BEAM LIMITING SYSTEMS

The following test procedure describes a method of checking simultaneously all the characteristics listed in 5.3.1.

- Step 1: Check that the PATIENT SUPPORT is horizontal using a spirit level. If a built-in angulation-indicator is used, check it against the spirit level and adjust if necessary. If the built-in angulation-indicator cannot be adjusted, record the deviation.*
- If the PATIENT SUPPORT cannot be brought to the horizontal position, record the deviation.*
- Step 2: Align the X-RAY SOURCE ASSEMBLY with the centre of the PATIENT SUPPORT or IMAGE RECEPTION AREA and adjust the X-RAY TUBE ASSEMBLY axis to be parallel to the long axis of the PATIENT SUPPORT using the facilities and indicators available on the equipment. Select the same FOCAL SPOT TO IMAGE RECEPTOR DISTANCE as used in the initial CONSTANCY TEST.*
- Step 3: Place a loaded 24 cm × 30 cm RADIOGRAPHIC CASSETTE (cassette "P": see figures 5a, 5b and 5c) in the centre of the CASSETTE CHANGER as determined by the clamps. Align the centre of the CASSETTE CHANGER with the X-RAY TUBE ASSEMBLY by means of the facilities provided and lock into position. Insert the CASSETTE CHANGER into the correct position for the X-RAY IMAGE RECEPTOR.*
- Step 4: For a manually adjusted or fixed BEAM LIMITING SYSTEM proceed directly to step 5.*
- Support a loaded 24 cm × 30 cm RADIOGRAPHIC CASSETTE (cassette "Q": see figure 5b) in a plane parallel to the PATIENT SUPPORT or X-RAY IMAGE RECEPTOR, but at a distance of between 20 cm and 30 cm closer to the X-RAY TUBE ASSEMBLY than the PATIENT SUPPORT or X-RAY IMAGE RECEPTOR. Orient the cassette so that the normal ENTRANCE SURFACE is on the side opposite the X-RAY TUBE ASSEMBLY and ensure that the cassette intercepts the entire X-RAY BEAM.*
- Measure and record the distance from the FOCAL SPOT to cassette "Q" at the centre of the LIGHT FIELD by means of the LIGHT FIELD INDICATOR and align the edges of the cassette parallel with the edges of the LIGHT FIELD.*

**NOTE** – The distances from the FOCAL SPOT to the two RADIOGRAPHIC CASSETTES are different and there are more absorbing materials in front of cassette "P", such as the PATIENT SUPPORT and the ANTI-SCATTER GRID and the cassette "Q". The AIR KERMA to the two cassettes may therefore differ considerably and, in an attempt to facilitate the choice of a single set of LOADING FACTORS to provide adequate IRRADIATION of the films in both cassettes, cassette "Q" is used in the reverse orientation to that normally employed.

If the test cannot satisfactorily be performed with a single IRRADIATION, it will be necessary to perform two IRRADIATIONS with different LOADING FACTORS.

- Step 5: For a manually adjusted or fixed BEAM LIMITING SYSTEM, place the alignment TEST DEVICE on the PATIENT SUPPORT or ENTRANCE SURFACE of the X-RAY IMAGE RECEPTOR. If the PATIENT SUPPORT is curved, place a thin, flat rigid surface of low ATTENUATION material, for example wood, across the PATIENT SUPPORT to maintain the TEST DEVICE in a plane parallel to the X-RAY IMAGE RECEPTOR.*

*For automatically adjusted BEAM LIMITING SYSTEMS place the alignment TEST DEVICE on the surface of cassette "Q" which is closest to the X-RAY TUBE ASSEMBLY.*

*The alignment TEST DEVICE incorporates the radio-opaque markings required in the evaluation of the size of the RADIATION FIELD.*

- Step 6:** *Place the alignment TEST DEVICE in the centre of the LIGHT FIELD according to the LIGHT FIELD INDICATOR and align the edges of the TEST DEVICE parallel to the edges of the LIGHT FIELD.*
- *For automatically adjusted BEAM LIMITING SYSTEMS record the size and position of the LIGHT FIELD in relation to the markings on the external surface of the TEST DEVICE.*
  - *For manually adjusted or fixed BEAM LIMITING SYSTEMS adjust the LIGHT FIELD size to the small field size markings on the alignment TEST DEVICE, for example 15 cm × 20 cm. Note the settings of the field size indicators of the BEAM LIMITING DEVICE. Record the extent of any asymmetry in the LIGHT FIELD, whereby one or more edges of the LIGHT FIELD cannot be set to coincide with the edges of the field defined by the alignment TEST DEVICE.*
- Step 7:** *Place the TEST DEVICE for the perpendicular position on top of the alignment TEST DEVICE so that the centre of the rings on the lower base of the TEST DEVICE for the perpendicular position coincides with the centre of the LIGHT FIELD indicated on the surface of the alignment TEST DEVICE. The test arrangements are shown in figures 5a and 5b.*
- Step 8:** *Make a radiographic IRRADIATION using the same LOADING FACTORS as in the initial CONSTANCY TEST in order to provide an optical density in the range from 0,5 to 1,5 on the processed RADIOGRAPHIC FILM.*

**NOTE** – *If the IRRADIATION can only be made under AUTOMATIC EXPOSURE CONTROL, it may also be necessary to use the ATTENUATION PHANTOM (see 5.1.2) as a substitute for the PATIENT, in order to provide appropriate ATTENUATION and hardening of the X-RAY BEAM interacting with the AUTOMATIC CONTROL SYSTEM. For the placement of the ATTENUATION PHANTOM, see figures 5a and 5b.*

- *For automatically adjusted BEAM LIMITING SYSTEMS, process the films in both cassettes.*
  - *For manually adjusted or fixed BEAM LIMITING SYSTEMS, re-adjust the LIGHT FIELD to the maximum cassette size. Process the film in the cassette.*
- Step 9:** *Remove the TEST DEVICE and repeat the above procedure (steps 1 to 8) but using loaded 35 cm × 43 cm RADIOGRAPHIC CASSETTES as cassettes "P" and "Q".*

**NOTE** – *A smaller RADIOGRAPHIC CASSETTE may be adequate at location "Q".*

*In addition, for fixed/manually adjusted BEAM LIMITING SYSTEMS, in step 6 use the larger field settings indicated on the alignment TEST DEVICE (for example 30 cm × 40 cm).*

### **5.3.3.3 X-RAY EQUIPMENT in which the RADIATION BEAM can only be directed horizontally**

*The procedures described in items 5.3.3.1 and 5.3.3.2 above are applicable to over-couch X-RAY TUBE ASSEMBLIES. The same procedure may be adopted for X-RAY EQUIPMENT in which the RADIATION BEAM is directed horizontally, except that the procedures will need to be modified as follows.*

- *Check that the PATIENT SUPPORT or X-RAY IMAGE RECEPTOR is in the vertical position instead of the horizontal position, using the built-in angulation-indicator or a spirit level if no such built-in angulation-indicator is provided.*
- *Provide means of supporting the TEST DEVICES in position versus the vertical surface of the PATIENT SURFACE or X-RAY IMAGE RECEPTOR input surface.*

#### 5.3.3.4 X-RAY EQUIPMENT with under-couch X-RAY TUBE ASSEMBLIES

*Because there is no light field-indication available with under-couch X-RAY TUBE ASSEMBLIES, it is not possible to perform checks on all the geometric characteristics described in 5.3.1. The only characteristic which can be checked without using methods of considerable complexity is the coincidence of the RADIATION FIELD with the X-RAY IMAGE RECEPTOR.*

**Step 1:** Check the PATIENT SUPPORT to be horizontal using the built-in angulation-indicator or a spirit level if no built-in angulation-indicator is provided. If the PATIENT SUPPORT is not horizontal or cannot be brought to the horizontal position, record the deviation.

**Step 2:** The under-couch X-RAY SOURCE ASSEMBLY will normally be aligned automatically with the centre of the IMAGE RECEPTION AREA when the X-RAY IMAGE RECEPTOR is in the operational position. Select the same FOCAL SPOT TO IMAGE RECEPTOR DISTANCE as used in the initial CONSTANCY TEST.

*For convenience this may be the maximum height at which the X-RAY IMAGE RECEPTOR can be set above the PATIENT SUPPORT.*

*Remove any secondary BEAM LIMITING DEVICE between the PATIENT and image receptor.*

**Step 3:** Place a loaded 24 cm × 30 cm RADIOGRAPHIC CASSETTE (cassette "P") in the centre of the CASSETTE CHANGER as determined by the clamps. Insert the CASSETTE CHANGER into the correct position for the X-RAY IMAGE RECEPTOR. The test arrangement is shown in figure 5c.

**Step 4:** Proceed to carry out this test in accordance with the method described for automatically adjusted BEAM LIMITING SYSTEMS in steps 4, 5, 8, and 9 of 5.3.3.2.

### 5.3.4 Data evaluation

#### 5.3.4.1 Manually adjusted or fixed BEAM LIMITING SYSTEMS

**Step 1:** Measure the distances between the edges of the RADIATION FIELD as defined by the darker irradiated area on each RADIOGRAM produced in cassette "P" and the edges of the LIGHT FIELD as indicated by the images of the radio-opaque markers. For each pair of opposite edges add these measurements to calculate the degree of misalignment; see figure 6a.

**Step 2:** Measure the lengths of the edges of the RADIATION FIELD as determined by the darker irradiated area on each of the RADIOGRAMS produced in cassette "P". Subtract the indicated dimensions of field size from the measured dimensions; see figure 6b.

#### 5.3.4.2 Automatic BEAM LIMITING SYSTEMS except those with under-couch X-RAY TUBE ASSEMBLIES

**Step 1:** Measure the distances between the edges of the RADIATION FIELD as defined by the darker irradiated area on each RADIOGRAM produced in cassette "Q" and the edges of the LIGHT FIELD as indicated by the images of the radio-opaque markers which correspond to the size and position of the LIGHT FIELD noted during the test; see step 2 of 5.3.3.2. Add these measurements for each pair of opposite edges to calculate the degree of misalignment.

**Step 2:** Measure the distances  $d_P$  and  $d_Q$  between two radio-opaque markers visible in both RADIOGRAMS, produced in cassettes "P" and "Q". Determine the magnification "m" by calculating the ratio  $d_P/d_Q$ .

*Measure the distances between the edges of the RADIATION FIELD as defined by the darker irradiated area on the RADIOGRAM produced in cassette "Q". To determine the distances between the edges of the RADIATION FIELD at the position of cassette "P", multiply these measured distances by "m".*

*Determine the difference between the estimated distance between the opposite edges of the RADIATION FIELD and the corresponding distance between the edges of the IMAGE RECEPTION AREA. Add these differences for each pair of opposite edges to calculate the degree of misalignment.*

### 5.3.4.3 Under-couch X-RAY TUBE ASSEMBLIES

*Make the measurements described in step 2 of 5.3.4.2.*

NOTE – This test cannot be performed if any edge of the RADIATION FIELD is outside the boundary of the X-RAY IMAGE RECEPTOR.

## 5.3.5 Criteria to be applied

### 5.3.5.1 Indicated FOCAL SPOT TO IMAGE RECEPTOR DISTANCE

The FOCAL SPOT TO IMAGE RECEPTOR DISTANCE shall be within  $\pm 1\%$  of the indicated distance and within  $\pm 1\%$  of the distance measured in the initial CONSTANCY TEST.

### 5.3.5.2 Perpendicular position of RADIATION BEAM AXIS to X-RAY IMAGE RECEPTOR (see note in 5.3.1)

The RADIATION BEAM AXIS shall be within  $1,5^\circ$  of the perpendicular axis to the IMAGE RECEPTION AREA. To check this, the image of the centre of the cross shall fall within the image of the inner circle of the TEST DEVICE for the perpendicular position.

### 5.3.5.3 Coincidence of RADIATION FIELD WITH LIGHT FIELD

In figure 6a, the measured discrepancies are represented by  $a_1$  and  $a_2$  on one axis and by  $b_1$  and  $b_2$  on the other. If the distance from the FOCAL SPOT is  $S$ , then, for compliance, the following relationships are true:

$$|a_1| + |a_2| \leq 0,02 \times S$$

$$|b_1| + |b_2| \leq 0,02 \times S$$

### 5.3.5.4 Coincidence of RADIATION FIELD WITH X-RAY IMAGE RECEPTOR

In figure 6b, the measured discrepancies are represented by  $c_1$  and  $c_2$  on one axis and by  $d_1$  and  $d_2$  on the other. If the distance from the FOCAL SPOT is  $S$ , then, for compliance, the following relationships are true:

$$|c_1| + |c_2| \leq 0,03 \times S$$

$$|d_1| + |d_2| \leq 0,03 \times S$$

$$|c_1| + |c_2| + |d_1| + |d_2| \leq 0,04 \times S$$

If, in addition to automatic beam limitation, a secondary BEAM LIMITING DEVICE is normally located between the PATIENT and the X-RAY IMAGE RECEPTOR, the misalignment criteria shall apply to the projection of the RADIATION FIELD as though the secondary BEAM LIMITING DEVICE were not present.

### 5.3.5.5 Accuracy of the numerical indication of the RADIATION FIELD size

The difference between the indicated dimensions of the RADIATION FIELD size and the measured dimensions shall be within  $\pm 2\%$  of the FOCAL SPOT TO IMAGE RECEPTOR DISTANCE.

### 5.3.6 Action to be taken

If the system fails to meet the criteria, the guidance given in annex C should be followed.

### 5.3.7 Frequency of CONSTANCY TESTS

The tests shall be carried out according to the INSTRUCTIONS FOR USE as provided by the MANUFACTURER. If no such information is provided, the CONSTANCY TESTS shall be performed at least quarterly.

## 5.4 Resolution of high-contrast detail

### 5.4.1 Summary

If variation in FOCAL SPOT size with time cannot be excluded, the following tests shall be performed.

This test checks the constancy of resolution of the X-RAY EQUIPMENT by producing a radiographic image of a high-contrast TEST DEVICE.

### 5.4.2 Test equipment

The following test equipment is required:

- a magnifying glass. A  $2,5 \times$  magnification is usually suitable;
- a high-contrast TEST DEVICE containing periodic patterns of radio-opaque materials (a detailed description of the high-contrast TEST DEVICE is provided in annex D);
- a dedicated RADIOGRAPHIC CASSETTE and associated screen-film combination of the type normally used with the X-RAY EQUIPMENT in clinical practice.

### 5.4.3 Test procedure

*Place the dedicated loaded cassette in the IMAGE RECEPTOR PLANE at the FOCAL SPOT TO IMAGE RECEPTOR DISTANCE used in normal clinical practice.*

*Place the high-contrast TEST DEVICE on the PATIENT SUPPORT, or the ENTRANCE PLANE of the SPOT FILM DEVICE so that the bar pattern is kept at about 200 mm from the surface (see annex D). Centre the X-RAY TUBE ASSEMBLY to the X-RAY IMAGE RECEPTOR and centre the test pattern to the X-RAY TUBE ASSEMBLY.*

*Orientate the main directions of the bars of the test pattern by approximately  $45^\circ$  to the X-RAY TUBE axis. For any item of test equipment, choose the measurement geometry so that the distance from the FOCAL SPOT and the position within the RADIATION FIELD can be reproduced to within  $\pm 1\%$  of the distance and position used in the initial CONSTANCY TEST.*

*Set the X-RAY TUBE VOLTAGE at about 70 kV and choose the LOADING FACTORS so as to provide an optical density of between 0,7 and 1,3 above FILM BASE PLUS FOG DENSITY in the unattenuated region of the processed film. Use the same LOADING FACTORS as those used in the initial CONSTANCY TEST.*

*NOTE – If the IRRADIATION can only be made under AUTOMATIC EXPOSURE CONTROL, it may also be necessary to use an ATTENUATION PHANTOM (see 5.1.2) as a substitute for the PATIENT, in order to provide appropriate ATTENUATION and hardening of the X-RAY BEAM interacting with the AUTOMATIC CONTROL SYSTEM. Under these circumstances, the ATTENUATION PHANTOM should be placed in the X-RAY BEAM as close as possible to the FOCAL SPOT as the X-RAY SOURCE ASSEMBLY will permit.*

*Take a RADIOGRAM with an X-RAY BEAM size of 100 mm  $\times$  100 mm measured at the PATIENT SUPPORT plane. Repeat the procedure for each available FOCAL SPOT size which can be selected.*

*For each test, note the LOADING FACTORS.*

*Process the irradiated films in accordance with the procedure referred to in 4.4.*

#### **5.4.4 Data evaluation**

*Examine the RADIOGRAM with the aid of a magnifying glass and record the maximum spatial frequencies visible. These are the cut-off frequencies under the test conditions.*

*In order to provide a comparison, undertake the same procedure with the RADIOGRAMS from the initial CONSTANCY TEST each time that this test is performed.*

#### **5.4.5 Criteria to be applied**

By comparison with the cut-off frequency from the initial CONSTANCY TEST, the measured cut-off frequency shall not decrease by more than:

- 20 % for continuously varying resolution test patterns; or
- one line pair group.

#### **5.4.6 Action to be taken**

If the system fails to meet the criteria, the guidance given in annex C should be followed.

#### **5.4.7 Frequency of CONSTANCY TESTS**

The tests shall be carried out according to the INSTRUCTIONS FOR USE as provided by the MANUFACTURER. If no such information is provided, the CONSTANCY TESTS shall be performed at least quarterly.

### **5.5 Variation in optical density throughout a RADIOGRAM**

#### **5.5.1 Summary**

This test checks the constancy of the relative variations in the optical density measured at specific points of a RADIOGRAM. This test may be performed simultaneously with the test described in 5.2, i.e. one single RADIOGRAM provides the information required for both tests.

#### **5.5.2 Test equipment**

Measurements are performed using the same dedicated RADIOGRAPHIC CASSETTE and associated screen-film combination of the type normally used with the X-RAY EQUIPMENT under test. Use a densitometer which reads consistently within  $\pm 0,02$ .

An ATTENUATION PHANTOM is used as a substitute for the PATIENT. This is intended to provide appropriate ATTENUATION and hardening of the X-RAY BEAM. In addition, a film marker TEST DEVICE shall be used in order to allow for the identification of specific points in the RADIOGRAM where the measurements of the optical density are to be made; see figure 1. Detailed descriptions of an ATTENUATION PHANTOM and a film marker TEST DEVICE are provided in annex D.

### 5.5.3 Test procedure

*Align the X-RAY SOURCE ASSEMBLY with the X-RAY IMAGE RECEPTOR as in normal clinical practice.*

*Place the loaded RADIOGRAPHIC CASSETTE in the IMAGE RECEPTION AREA and the film marker TEST DEVICE in the X-RAY BEAM between the FOCAL SPOT and the X-RAY IMAGE RECEPTOR.*

*If IRRADIATIONS can only be made under AUTOMATIC EXPOSURE CONTROL, use an ATTENUATION PHANTOM as a substitute for the PATIENT in order to provide appropriate ATTENUATION and hardening of the X-RAY BEAM.*

*For any item of test equipment, the measurement geometry shall be the same within  $\pm 1$  % as used in the initial CONSTANCY TEST. Use the same RADIATION FIELD size.*

*For the initial CONSTANCY TEST select at least two X-RAY TUBE VOLTAGE settings that are typical for the radiological procedures for which the X-RAY EQUIPMENT is used. For the following CONSTANCY TESTS select the same settings as those used in the initial CONSTANCY TEST.*

#### 5.5.3.1 Testing under manual control

*Operate the X-RAY EQUIPMENT using manual settings for the LOADING FACTORS identical to those used in the initial CONSTANCY TEST.*

*Process the irradiated film in accordance with the procedure referred to in 4.4.*

*Measure the optical density at the specific points on the RADIOGRAM indicated by the images of the radio-opaque markers; see figure 1.*

*If necessary, adjust the measured values of density in accordance with the procedure referred to in 4.4 to allow for any changes in film batch or processing conditions.*

#### 5.5.3.2 Testing under AUTOMATIC EXPOSURE CONTROL

*Operate the X-RAY EQUIPMENT in conjunction with the AUTOMATIC CONTROL SYSTEM and use the X-RAY TUBE VOLTAGE settings and all other appropriate settings, for example density control, identical to those used in the initial CONSTANCY TEST.*

*Process the irradiated film in accordance with the procedure referred to in 4.4.*

*Measure the optical density at the specific points on the RADIOGRAM indicated by the images of the radio-opaque markers; see figure 1.*

*If necessary, adjust the measured values of density in accordance with the procedure referred to in 4.4 to allow for any changes in film batch or processing conditions.*

### 5.5.4 Data evaluation

*Examine the RADIOGRAM and compare it with that produced in the initial CONSTANCY TEST. This will reveal gross variations in the distribution of optical density throughout the RADIOGRAM.*

*Determine the differences in optical density measured at specific points on the RADIOGRAM, indicated by the images of the radio-opaque markers, and that measured at the reference point, and compare with the established initial CONSTANCY TEST. Make additional measurements and determine the differences in optical density at any other points on the RADIOGRAM where visual inspection indicates that a significant change may have occurred.*

### 5.5.5 Criteria to be applied

The optical density differences should be within  $\pm 0,10$  of the BASELINE VALUES.

### 5.5.6 Action to be taken

If the system fails to meet the criteria, the guidance given in annex C should be followed.

### 5.5.7 Frequency of variation in optical densities CONSTANCY TESTS

The tests shall be carried out according to the INSTRUCTIONS FOR USE as provided by the MANUFACTURER. If no such information is provided, the CONSTANCY TESTS shall be performed at least quarterly.

## 6 Statement of compliance

The test report shall be headed


**Test report  
on constancy test of equipment for general direct radiography  
according to IEC 61223-2-11:1999**

If compliance with this standard is to be stated, this shall be done as follows:


The equipment for general direct radiography, (...\*), complies with IEC 61223-2-11:1999.

---

\* Identification (for example, name of equipment, model or type reference).

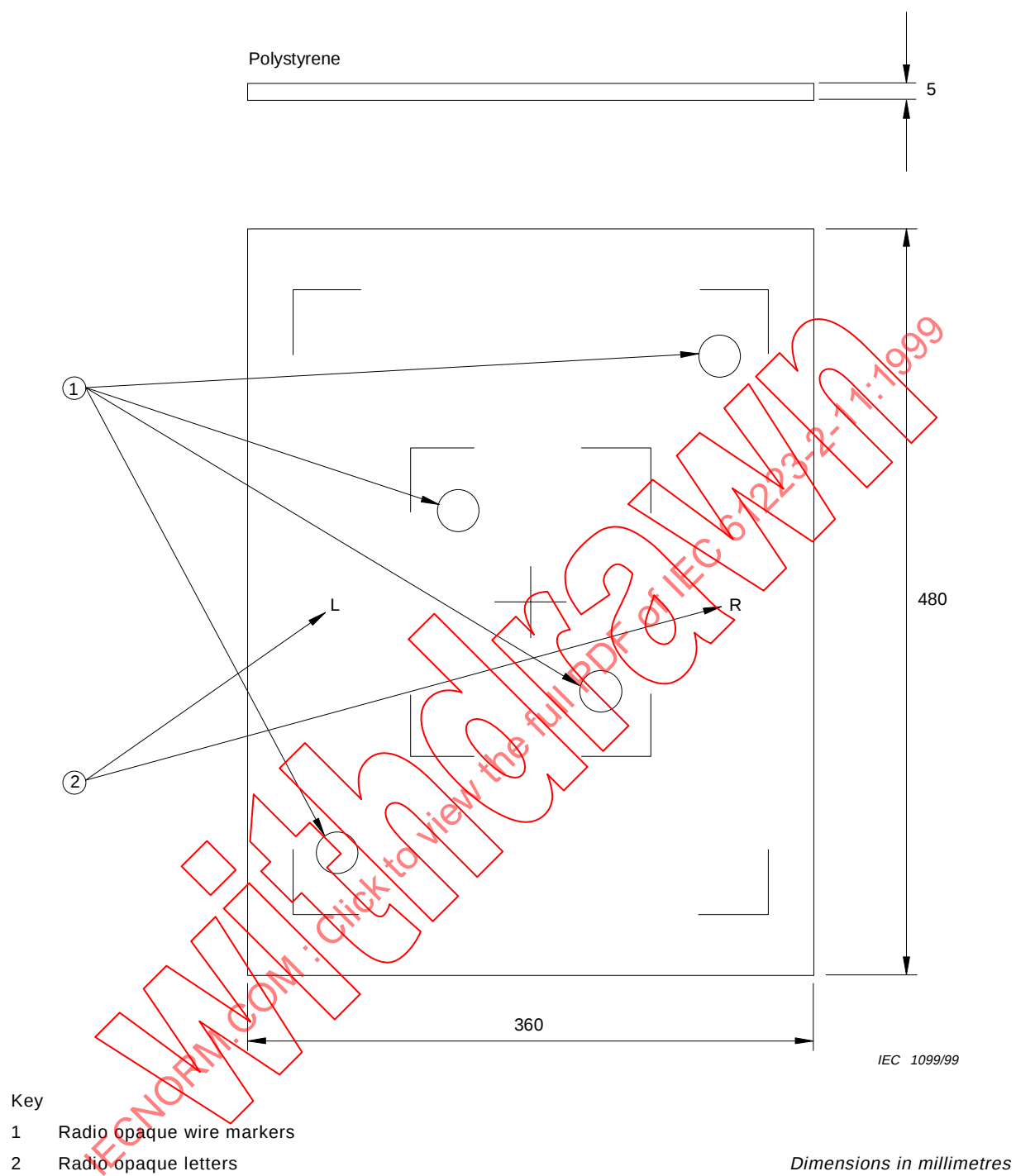

**Figure 1 – Example of a film marker TEST DEVICE**

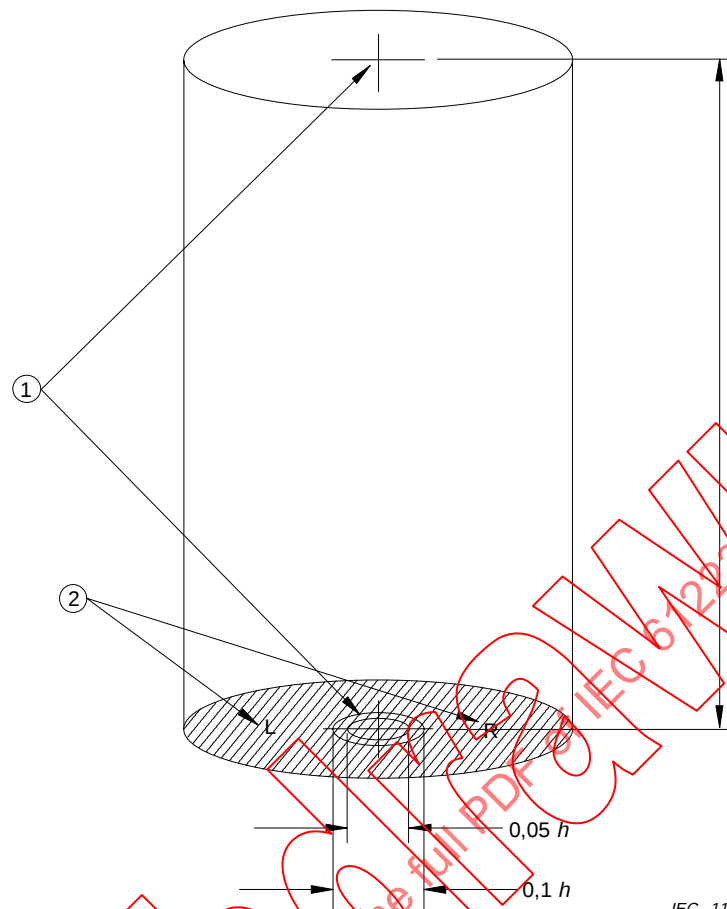


IEC 1100/99

Key

- 1 Radio opaque wire markers
- 2 Radio opaque letters

*Dimensions in millimetres*

**Figure 2 – Example of a TEST DEVICE for perpendicular position**

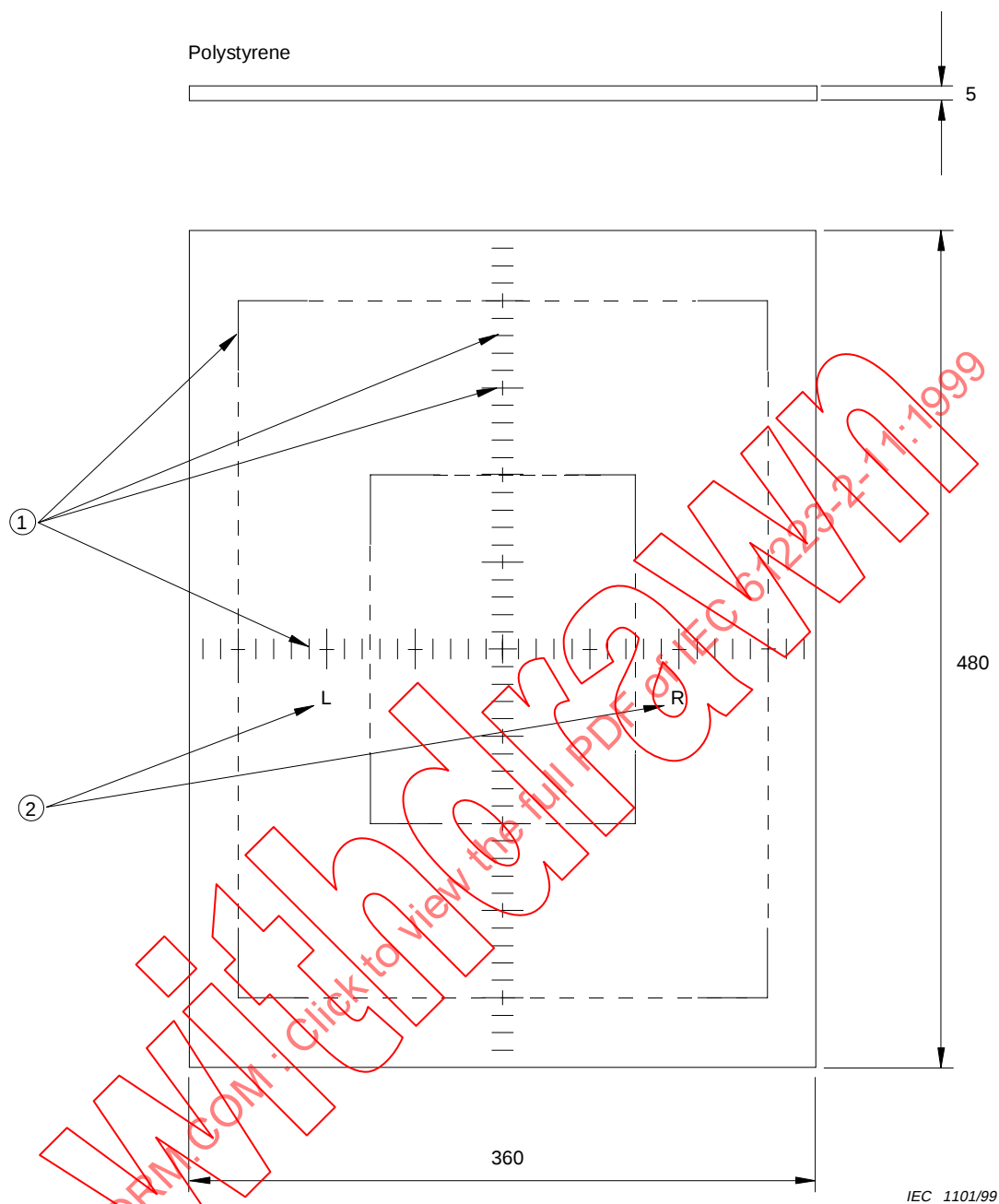

## Key

- 1 Radio opaque wire markers
- 2 Radio opaque letters

*Dimensions in millimetres*

**Figure 3 – Example of an alignment TEST DEVICE**


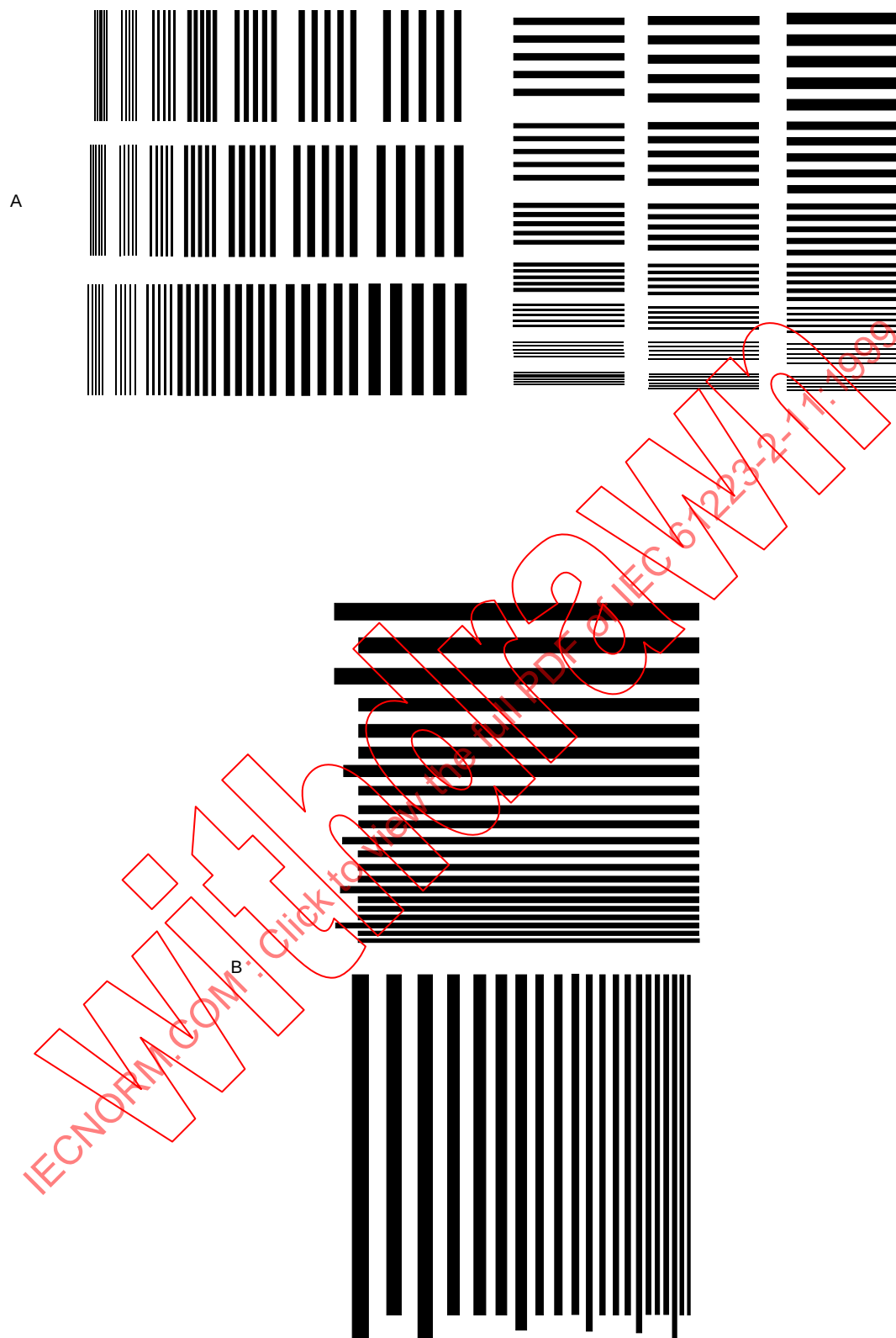


Figure 4 – Examples of a high-contrast TEST DEVICE

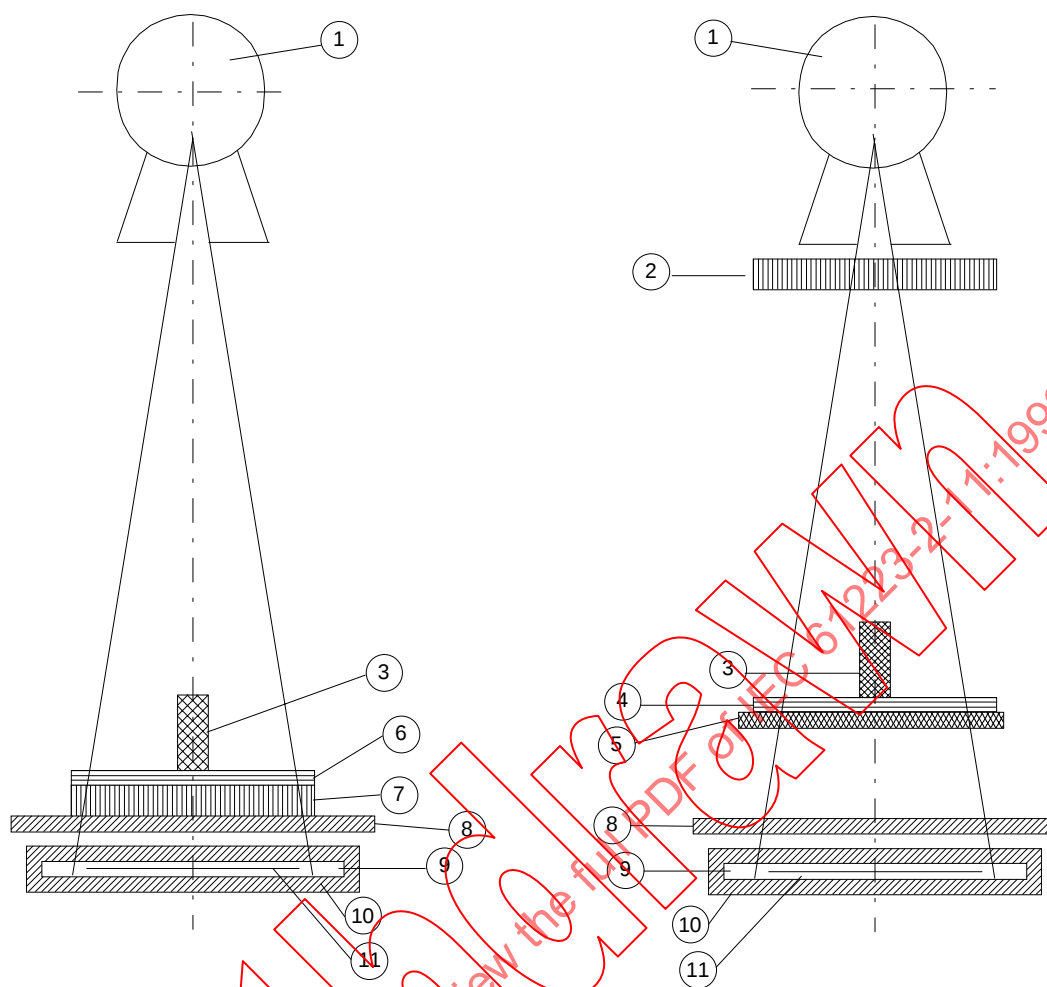


IEC 1103/99

IEC 1104/99

**Figure 5a – Over-couch arrangement****Figure 5b – Over-couch arrangement****Key**

- |                                          |                         |
|------------------------------------------|-------------------------|
| 1 X-RAY SOURCE ASSEMBLY                  | 7 ATTENUATION PHANTOM   |
| 2 ATTENUATION PHANTOM                    | 8 PATIENT SUPPORT       |
| 3 TEST DEVICE for perpendicular position | 9 CASSETTE P            |
| 4 Alignment TEST DEVICE                  | 10 X-RAY IMAGE RECEPTOR |
| 5 CASSETTE Q                             | 11 IMAGE RECEPTOR AREA  |
| 6 Alignment TEST DEVICE                  |                         |

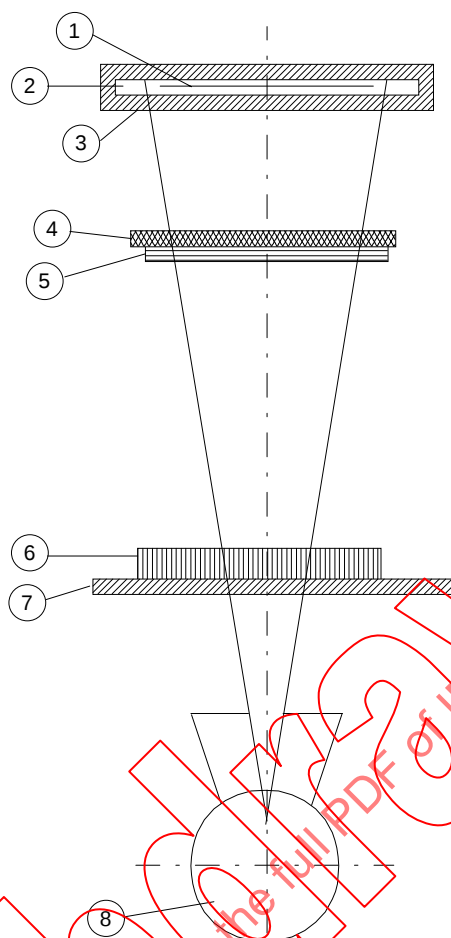



IEC 1105/99

**Figure 5c – Under-couch arrangement**

**Key**

- |   |                      |   |                       |
|---|----------------------|---|-----------------------|
| 1 | IMAGE RECEPTOR AREA  | 5 | ALIGNMENT TEST DEVICE |
| 2 | CASSETTE P           | 6 | ATTENUATION PHANTOM   |
| 3 | X-RAY IMAGE RECEPTOR | 7 | PATIENT SUPPORT       |
| 4 | CASSETTE Q           | 8 | X-RAY SOURCE ASSEMBLY |

**Figure 5 – Arrangements for testing geometric conditions**

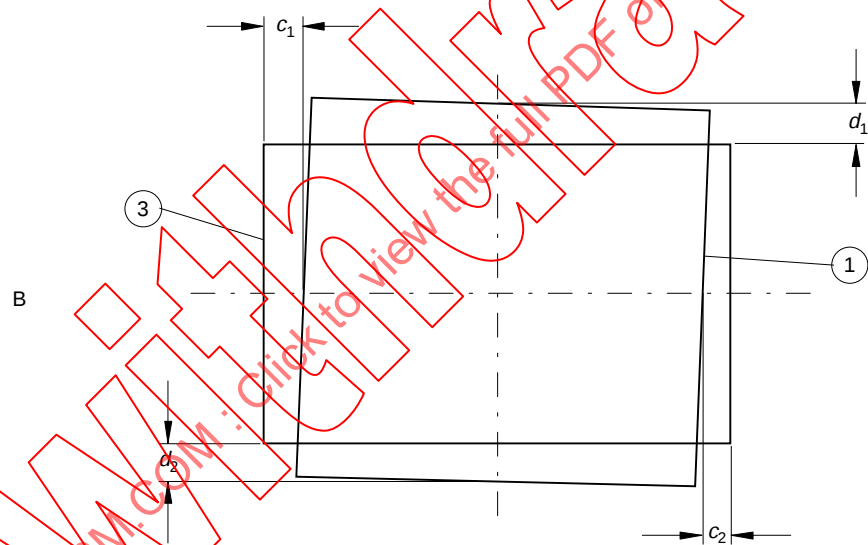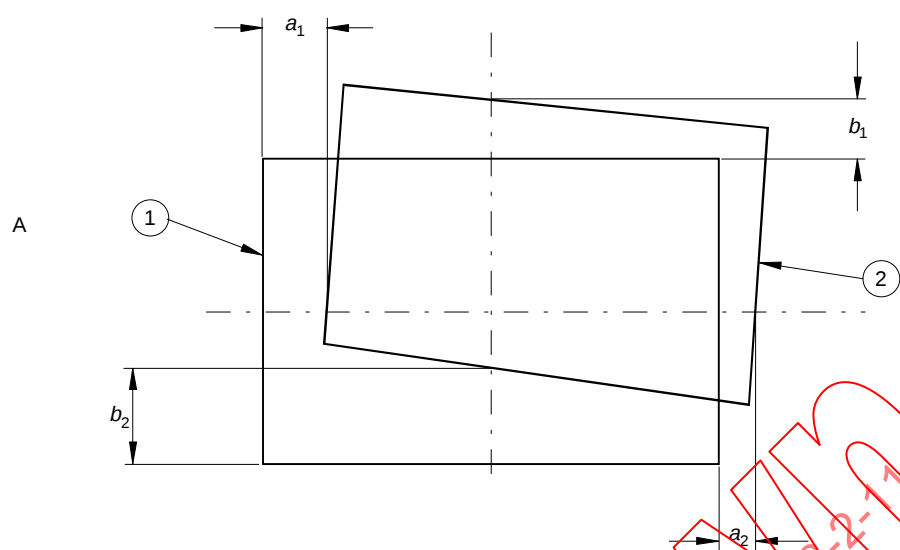

IEC 1106/99

## Key

- 1 RADIATION FIELD
- 2 LIGHT FIELD
- 3 IMAGE RECEPTION AREA

Figure 6 – Geometric coincidences

## Annex A (normative)

### Terminology – Index of defined terms

|                                                   |           |
|---------------------------------------------------|-----------|
| IEC 60788 .....                                   | rm.-.-..  |
| Name of unit in the International System SI ..... | rm.-.-.*  |
| Derived term without definition .....             | rm.-.-.+  |
| Term without definition .....                     | rm.-.-.-  |
| Name of earlier unit .....                        | rm.-.-.   |
| Shortened term .....                              | rm.-.-.s  |
| Clause 3 of IEC 61223-1 .....                     | AG-3...   |
| Clause 3 of IEC 61223-2-XY .....                  | XY-3...   |
| ACCEPTANCE TEST .....                             | AG-3.2.4  |
| ACCESSORY .....                                   | rm-83-06  |
| ACCOMPANYING DOCUMENTS .....                      | rm-82-01  |
| ADDED FILTER .....                                | rm-35-02  |
| AIR KERMA .....                                   | rm-13-11  |
| ANODE .....                                       | rm-22-06  |
| ANTI-SCATTER GRID .....                           | rm-32-06  |
| ATTENUATION .....                                 | rm-12-08  |
| AUTOMATIC CONTROL SYSTEM .....                    | rm-36-45  |
| AUTOMATIC EXPOSURE CONTROL .....                  | rm-36-46  |
| BASELINE VALUE .....                              | AG-3.2.7  |
| BEAM LIMITING DEVICE .....                        | rm-37-28  |
| BEAM LIMITING SYSTEM .....                        | rm-37-27  |
| CASSETTE CHANGER .....                            | rm-31-06  |
| CONSTANCY TEST .....                              | AG-3.2.6  |
| CURRENT TIME PRODUCT .....                        | rm-36-13  |
| DIRECT RADIOGRAPHY .....                          | rm-41-07  |
| ENTRANCE PLANE .....                              | rm-32-42  |
| ENTRANCE SURFACE .....                            | rm-37-17  |
| ESTABLISHED CRITERION .....                       | AG-3.2.8  |
| FILM BASE PLUS FOG DENSITY .....                  | 1-3.2.2   |
| FILM ILLUMINATOR .....                            | 2-3.2.1   |
| FILM PROCESSOR .....                              | 1-3.2.1   |
| FILTRATION .....                                  | rm-12-11  |
| FOCAL SPOT .....                                  | rm-20-13s |
| FOCAL SPOT TO IMAGE RECEPTOR DISTANCE .....       | rm-37-13  |
| HIGH-VOLTAGE GENERATOR .....                      | rm-21-01  |
| IMAGE RECEPTION AREA .....                        | rm-37-16  |
| IMAGE RECEPTOR PLANE .....                        | rm-37-15  |
| INSTRUCTIONS FOR USE .....                        | rm-82-02  |
| INTENSIFYING SCREEN .....                         | rm-32-38  |
| IONIZATION CHAMBER .....                          | rm-51-03  |
| IRRADIATION .....                                 | rm-12-09  |
| IRRADIATION TIME .....                            | rm-36-11  |

|                                   |           |
|-----------------------------------|-----------|
| LIGHT FIELD .....                 | rm-37-09  |
| LIGHT FIELD INDICATOR .....       | rm-37-31  |
| LOADING FACTOR .....              | rm-36-01  |
| LOADING TIME .....                | rm-36-10  |
| MANUFACTURER .....                | rm-85-03- |
| MOVING GRID .....                 | rm-32-15  |
| OPERATOR .....                    | rm-85-02  |
| PATIENT .....                     | rm-62-03  |
| PATIENT SURFACE .....             | rm-37-18  |
| PHANTOM .....                     | rm-54-01  |
| PRIMARY RADIATION .....           | rm-11-06  |
| PATIENT SUPPORT .....             | rm-30-02  |
| QUALITY ASSURANCE PROGRAMME ..... | AG-3.2.2  |
| RADIATION .....                   | rm-11-01  |
| RADIATION BEAM .....              | rm-37-05  |
| RADIATION BEAM AXIS .....         | rm-37-06  |
| RADIATION DETECTOR .....          | rm-51-01  |
| RADIATION FIELD .....             | rm-37-07  |
| RADIATION METER .....             | rm-50-01  |
| RADIOGRAM .....                   | rm-32-02  |
| RADIOGRAPHIC CASSETTE .....       | rm-35-14  |
| RADIOGRAPHIC FILM .....           | rm-32-32  |
| RADIOLOGICAL INSTALLATION .....   | rm-20-24  |
| REFERENCE DIRECTION .....         | rm-37-02  |
| SCATTERED RADIATION .....         | rm-11-13  |
| Specific .....                    | rm-74-01  |
| Specified .....                   | rm-74-02  |
| SPOTFILM DEVICE .....             | rm-31-05  |
| STATUS TEST .....                 | AG-3.2.5  |
| TARGET .....                      | rm-20-08  |
| TEST DEVICE .....                 | rm-71-04  |
| TISSUE EQUIVALENT MATERIAL .....  | rm-35-16  |
| TOTAL FILTRATION .....            | rm-13-48  |
| USER .....                        | rm-85-01  |
| X-RADIATION .....                 | rm-11-01- |
| X-RAY BEAM .....                  | rm-37-05+ |
| X-RAY EQUIPMENT .....             | rm-20-20  |
| X-RAY IMAGE RECEPTOR .....        | rm-32-29  |
| X-RAY SOURCE ASSEMBLY .....       | rm-20-05+ |
| X-RAY SPECTRUM .....              | rm-13-34+ |
| X-RAY TUBE .....                  | rm-22-03  |
| X-RAY TUBE ASSEMBLY .....         | rm-22-01  |
| X-RAY TUBE CURRENT .....          | rm-36-07  |
| X-RAY TUBE VOLTAGE .....          | rm-36-02  |