
**Plastics — Differential scanning
calorimetry (DSC) —**

**Part 4:
Determination of specific heat
capacity**

*Plastiques — Analyse calorimétrique différentielle (DSC) —
Partie 4: Détermination de la capacité thermique massique*



STANDARDSISO.COM : Click to view the full PDF of ISO 11357-4:2021



COPYRIGHT PROTECTED DOCUMENT

© ISO 2021

All rights reserved. Unless otherwise specified, or required in the context of its implementation, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
CP 401 • Ch. de Blandonnet 8
CH-1214 Vernier, Geneva
Phone: +41 22 749 01 11
Email: copyright@iso.org
Website: www.iso.org

Published in Switzerland

Contents

	Page
Foreword	iv
1 Scope	1
2 Normative references	1
3 Terms and definitions	1
4 Principle	2
4.1 General	2
4.2 Continuous-scanning method	3
4.3 Stepwise-scanning method	4
5 Apparatus	4
6 Test specimen	4
7 Test conditions and specimen conditioning	5
8 Procedure	5
8.1 Selection of crucibles	5
8.2 Setting up the apparatus and adjustment of isothermal baselines	5
8.3 Measurement of specific heat capacity of calibration material	6
8.4 Specimen run	8
9 Determination of specific heat capacities	8
9.1 Calculation of specific heat capacities	8
9.2 Numerical rounding of the results	8
10 Precision and bias	8
11 Test report	9
Annex A (informative) Approximate expression of the specific heat capacity of pure α-alumina [3][4]	10
Bibliography	12

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular, the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see www.iso.org/directives).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see www.iso.org/patents).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation of the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT), see www.iso.org/iso/foreword.html.

This document was prepared by Technical Committee ISO/TC 61, *Plastics*, Subcommittee SC 5, *Physical-chemical properties*, in collaboration with the European Committee for Standardization (CEN) Technical Committee CEN/TC 249, *Plastics*, in accordance with the Agreement on technical cooperation between ISO and CEN (Vienna Agreement).

This third edition cancels and replaces the second edition (ISO 11357-4:2014), which has been technically revised. The main changes compared to the previous edition are as follows:

- the measurement procedure has been updated;
- reference data of α -alumina have been updated.

A list of all parts in the ISO 11357 series can be found on the ISO website.

Any feedback or questions on this document should be directed to the user's national standards body. A complete listing of these bodies can be found at www.iso.org/members.html.

Plastics — Differential scanning calorimetry (DSC) —

Part 4:

Determination of specific heat capacity

1 Scope

This document specifies methods for determining the specific heat capacity of plastics by differential scanning calorimetry.

2 Normative references

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 472, *Plastics — Vocabulary*

ISO 11357-1, *Plastics — Differential scanning calorimetry (DSC) — Part 1: General principles*

ISO 80000-1, *Quantities and units — Part 1: General*

3 Terms and definitions

For the purposes of this document, the terms and definitions given in ISO 472, ISO 11357-1 and the following apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- ISO Online browsing platform: available at <https://www.iso.org/obp>
- IEC Electropedia: available at <http://www.electropedia.org/>

3.1

calibration material

material of known specific heat capacity

Note 1 to entry: Usually, α -alumina (such as synthetic sapphire) of 99,9 % or higher purity is used as the calibration material.

3.2

specific heat capacity at constant pressure

c_p

quantity of heat necessary to raise the temperature of unit mass of material by 1 K at constant pressure

Note 1 to entry: It is given by the following formula:

$$c_p = m^{-1} \cdot C_p = m^{-1} \cdot \left(\frac{dQ}{dT} \right)_p$$

where

- c_p is the specific heat capacity and is expressed in kilojoules per kilogram per K ($\text{kJ}\cdot\text{kg}^{-1}\cdot\text{K}^{-1}$) or in joules per gram per K ($\text{J}\cdot\text{g}^{-1}\cdot\text{K}^{-1}$); subscript p indicates an isobaric process;
- m is the mass of material, expressed in kilogram (kg) or gram (g);
- C_p is the total heat capacity and is expressed in kilojoules per K ($\text{kJ}\cdot\text{K}^{-1}$) or in joules per K ($\text{J}\cdot\text{K}^{-1}$); subscript p indicates an isobaric process;
- $\left(\frac{dQ}{dT}\right)_p$ is the quantity of heat dQ necessary to raise the temperature of the material by dT , expressed in kilojoules per K ($\text{kJ}\cdot\text{K}^{-1}$) or in joules per K ($\text{J}\cdot\text{K}^{-1}$), measured at constant pressure.

This formula is valid in a temperature range where a material shows no first-order phase transition.

The quotient $\left(\frac{dQ}{dT}\right)$ can be obtained by dividing the heat flow rate by the heating rate:

$$(dQ/dT) = \frac{(dQ/dt)}{(dT/dt)}$$

where

(dQ/dt) is the heat flow rate, expressed in kilojoules per second ($\text{kJ}\cdot\text{s}^{-1}$) or in joules per second ($\text{J}\cdot\text{s}^{-1}$) or in watts (W);

(dT/dt) is the heating rate, expressed in kelvins (K) per second (s) ($\text{K}\cdot\text{s}^{-1}$).

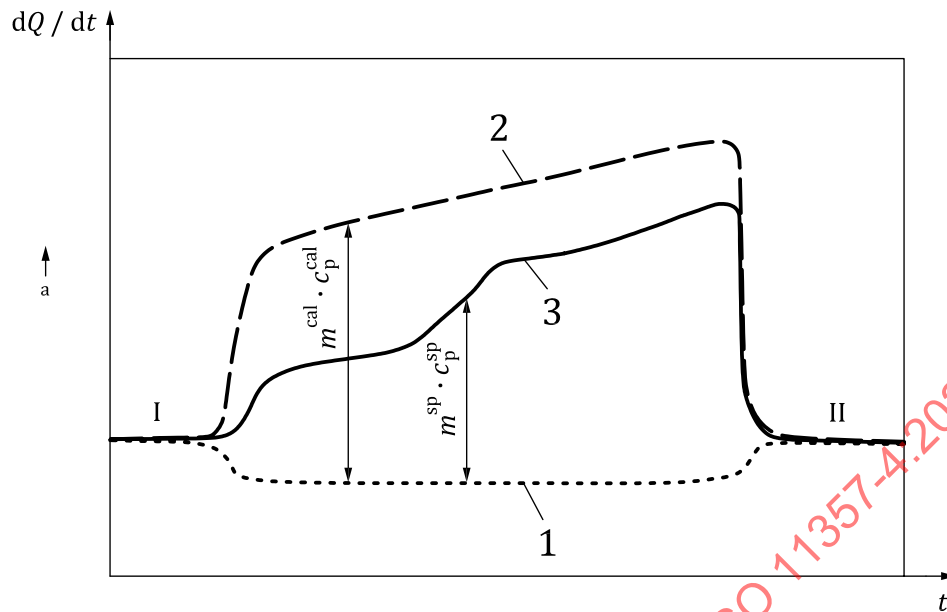
Note 2 to entry: At phase transitions, there is a discontinuity in the heat capacity. Part of the heat is consumed to produce a material state of higher energy and it is not all used in raising the temperature. For this reason, the specific heat can only be determined properly outside regions of phase transitions.

4 Principle

4.1 General

Each measurement consists of three runs at the same scanning rate (see [Figure 1](#)):

- a blank run (empty crucibles in sample and reference holders);
- a calibration run (calibration material in sample holder crucible and empty crucible in reference holder);
- a specimen run (specimen in sample holder crucible and empty crucible in reference holder).

**Key** dQ/dt heat flow rate t time

1 blank run

2 calibration run

3 specimen run

I isothermal baseline at start temperature, T_s II isothermal baseline at end temperature T_f

a Endothermic direction.

Figure 1 — Schematic drawing of typical DSC curves for specific heat capacity measurement (blank, calibration and specimen runs) after baseline adjustment

4.2 Continuous-scanning method

Based on the DSC principle (see ISO 11357-1) and the definition of specific heat capacity at constant pressure given in 3.2, the following relations shown as [Formulae \(1\)](#) and [\(2\)](#) can be obtained:

$$m^{\text{sp}} \cdot c_p^{\text{sp}} \propto P_{\text{sp}} - P_{\text{bl}} \quad (1)$$

$$m^{\text{cal}} \cdot c_p^{\text{cal}} \propto P_{\text{cal}} - P_{\text{bl}} \quad (2)$$

where

P is the heat flow rate (dQ/dt);

$^{\text{sp}}$ is the specimen;

$^{\text{bl}}$ is an empty crucible (blank run);

$^{\text{cal}}$ is the calibration material.

See [Figure 1](#).

When P_{sp} , P_{cal} and P_{bl} are measured, c_p^{sp} can be calculated using [Formula \(4\)](#), since the values of c_p^{cal} , m^{sp} and m^{cal} are known:

$$\frac{m^{sp} \cdot c_p^{sp}}{m^{cal} \cdot c_p^{cal}} = \frac{P_{sp} - P_{bl}}{P_{cal} - P_{bl}} \quad (3)$$

$$c_p^{sp} = c_p^{cal} \cdot \frac{m^{cal} (P_{sp} - P_{bl})}{m^{sp} (P_{cal} - P_{bl})} \quad (4)$$

4.3 Stepwise-scanning method

In the stepwise-scanning method, the total temperature range to be scanned is divided into small intervals and a complete determination consisting of the three runs mentioned in [4.1](#) is performed for each temperature interval. Upon integration of the heat flow rate curve, the total heat ΔQ consumed in the interval can be obtained. Dividing ΔQ by the temperature interval ΔT and the mass of the specimen gives the specific heat (see [3.2](#)):

$$m^{sp} \cdot c_p^{sp} \propto \left(\frac{\Delta Q^{sp}}{\Delta T} \right)_p - \left(\frac{\Delta Q^{bl}}{\Delta T} \right)_p \quad (5)$$

$$m^{cal} \cdot c_p^{cal} \propto \left(\frac{\Delta Q^{cal}}{\Delta T} \right)_p - \left(\frac{\Delta Q^{bl}}{\Delta T} \right)_p \quad (6)$$

Keeping the temperature intervals ΔT constant, combining [Formula \(5\)](#) and [Formula \(6\)](#) results in:

$$c_p^{sp} = c_p^{cal} \cdot \frac{m^{cal}}{m^{sp}} \cdot \frac{\Delta Q^{sp} - \Delta Q^{bl}}{\Delta Q^{cal} - \Delta Q^{bl}} \quad (7)$$

5 Apparatus

5.1 DSC apparatus. See ISO 11357-1.

5.2 Crucibles. See ISO 11357-1.

The crucibles for the test specimen and the reference specimen (calibration material) shall be of the same shape and material and their masses shall not differ by more than 0,1 mg.

NOTE The same blank run and calibration run can be used for several measurements, if the instrument is sufficiently stable and the difference in mass between the calibration material crucible and the empty crucible is corrected for. An adequate correction can be obtained by adding the term $c_{p,crucible}(T)\beta\Delta m$ to the heat flow rate of the calibration run, where $c_{p,crucible}(T)$ is the specific heat capacity of the calibration crucible as a function of temperature, β is the heating rate and Δm is the difference in mass between the calibration crucible and the empty crucible. The same procedure can also be used for correcting differences in mass between the specimen run and the blank run.

5.3 Analytical balance. See ISO 11357-1.

6 Test specimen

See ISO 11357-1.

7 Test conditions and specimen conditioning

See ISO 11357-1.

8 Procedure

8.1 Selection of crucibles

Prepare three crucibles and their lids and weigh the crucibles together with their lids. The total mass shall not differ by more than 0,1 mg (see 5.2). In other respects, such as material, size, crucible type (open or sealed), the crucibles shall be identical.

NOTE It is also possible to use the same crucible in the sample holder for blank run, calibration run, and specimen run if the lid is just laid upon the crucible for blank run and calibration run. For the specimen run, the crucible can be hermetically sealed.

8.2 Setting up the apparatus and adjustment of isothermal baselines

8.2.1 Place a pair of empty crucibles with lids in the DSC sample and reference holders.

8.2.2 If using a continuous-scan programme:

- a) Set the start and end temperatures (T_s and T_f). The start temperature T_s should be at least 30 K lower than that at which data are first required.

When more precise results are required over a wide temperature range, the overall range may be divided into two (or more) smaller ranges, each 50 K to 100 K wide. The start temperature T_s of the second range should be 30 K lower than the end temperature T_f of the first temperature range to ensure sufficient overlap.

- b) Set the scanning rate.
- c) Set the dwell time of the isothermal stages I and II (see Figure 1) and allow the respective isothermal baselines to stabilize. This dwell time will usually be between 2 min and 10 min.

NOTE Some calorimeters, such as those of the Calvet type, can require up to 30 min before the baseline stabilizes.

8.2.3 If using a stepwise-scan programme:

When the specific heat capacities of the samples do not significantly depend on the temperature, the stepwise-scanning method may be used in which the integration of the heat flow over small temperature intervals gives a set of individual specific heat values for the temperature intervals considered. Attention shall be paid to the following points.

- a) The dwell time of the isothermal stages shall be sufficiently long to obtain a stable baseline.
- b) This method shall not be used over a temperature range in which first-order phase transitions occur.

The stepwise scan is performed as follows.

- Set the start and end temperatures (T_s and T_f).
- Set the temperature increment preferably to 5 K or 10 K.
- Set the temperature-scanning rate to 5 K·min⁻¹ or 10 K·min⁻¹.
- Set the dwell time of the isothermal stages, usually to between 2 min and 10 min.

8.2.4 Set the sensitivity of the heat flow rate in order to obtain an ordinate span of at least 80 % of full scale (see [Figure 1](#)).

8.2.5 Adjust the apparatus so that the isothermal baselines before and after the heating stage are at the same ordinate level.

Check that adjustment of the baselines of the respective DSC curves results in the same ordinate level. If the baseline reproducibility is poor, readjust the apparatus and repeat the determination.

NOTE Other reasons for poor baseline reproducibility can be contamination of the sample crucible, the position of the lid, the stability of the purge gas flow rate, sample decomposition, sample evaporation, chemical reaction between crucible and sample, etc.

8.2.6 Execute the temperature programme set as described in [8.2.2](#) or [8.2.3](#). [Figure 2](#) shows a typical DSC curve obtained in the continuous-scanning mode whereas [Figure 3](#) shows a DSC curve obtained in the stepwise-scanning mode.

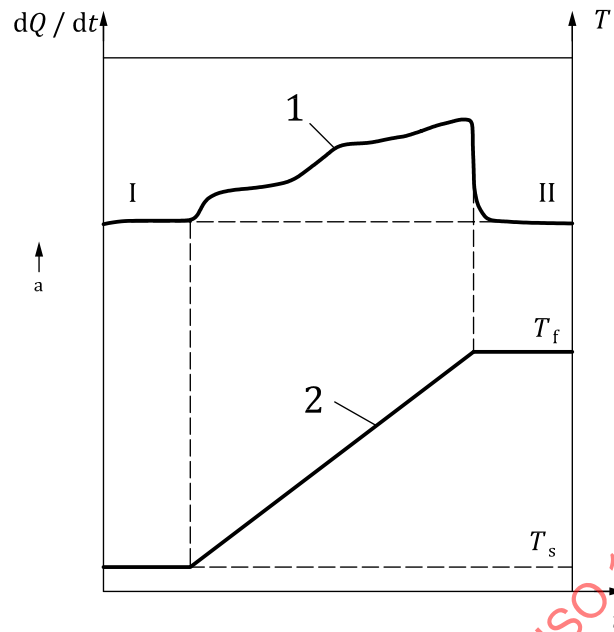
8.3 Measurement of specific heat capacity of calibration material

Using an analytical balance, weigh a calibration material, such as α -alumina (synthetic sapphire) of 99,9 % or higher purity, into one of the crucibles prepared in [8.1](#). Put the crucible containing the calibration material, with the lid, in the sample holder and carry out the same measurement(s) as described for the blank run(s) in [8.2](#).

NOTE Small differences in the masses of the crucibles used for the specimen, calibration and blank runs can be corrected for as indicated in the Note to [5.2](#).

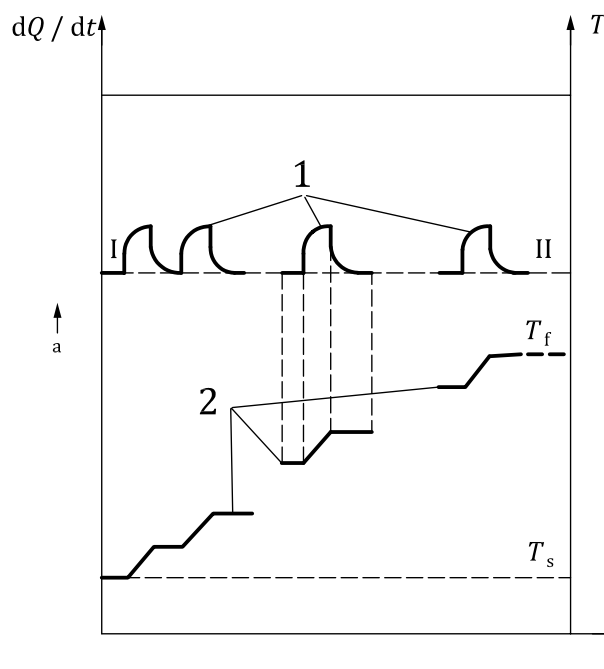
The heat capacity of the calibration material should match that of the specimen to be analysed as closely as possible in order to minimize systematic errors.

The nominal values of the specific heat capacity of α -alumina at various temperatures are given in [Annex A, Table A.1](#).

**Key**

- 1 DSC curve
- 2 temperature curve
- I isothermal baseline at start temperature T_s
- II isothermal baseline at end temperature T_f
- a Endothermic direction.

Figure 2 — Schematic drawing of a continuous-scan DSC curve

**Key**

- 1 DSC curve
- 2 temperature curve
- I isothermal baseline at start temperature T_s
- II isothermal baseline at end temperature T_f
- a Endothermic direction.

Figure 3 — Schematic drawing of a stepwise-scan DSC curve

8.4 Specimen run

Weigh the test specimen into the sample crucible. Put the crucible containing the test specimen, with the lid, into the sample holder and perform a DSC run as for the calibration material. A large mass of sample is recommended.

The same blank run as was used for the calibration run in 8.3 may also be used for the specimen run.

9 Determination of specific heat capacities

9.1 Calculation of specific heat capacities

Calculate c_p^{sp} , in $J \cdot g^{-1} \cdot K^{-1}$, using [Formula \(4\)](#) for continuous heating or [Formula \(7\)](#) for the stepwise method.

9.2 Numerical rounding of the results

Round the specific heat capacity values thus obtained to the second decimal place, using the method specified in ISO 80000-1.

10 Precision and bias

The precision and bias of the methods described are not known at the time of publication.

11 Test report

The following items shall be reported:

- a) a reference to this document, i.e. ISO 11357-4:2021;
- b) the date of the test;
- c) all details necessary for complete identification of the sample tested, including the thermal history;
- d) the manufacturer, model and type (power compensation or heat flux) of DSC apparatus used;
- e) the shape, dimensions and materials of the crucible and lid;
- f) the test atmosphere and the flow rate of the inflow gas;
- g) the calibration material, including information on the issuing body, the nature of the material, the mass used and other characteristics relevant to the calibration;
- h) the shape, dimensions and mass of the test specimen;
- i) details of sampling and the conditioning of the test specimen;
- j) the temperature programme parameters, i.e. the start temperature, the heating rate, the end temperature, the dwell time of the isothermal stages and, in the stepwise-heating method, the temperature increment, as well as the cooling rate in conditioning, if carried out;
- k) the test results, including the specific heat capacities and the respective temperatures;
- l) any deviations from the specified procedure;
- m) any unusual features observed;
- n) other items required, if any.

Annex A (informative)

Approximate expression of the specific heat capacity of pure α -alumina [3][4]

The specific heat capacity values of corundum (α -Al₂O₃, "Sapphire") in [Table A.2](#) are approximated by [Formula \(A.1\)](#):

$$c_p = a_0 + a_1 \cdot T + \frac{a_2}{T^2} + \sum_{n=3}^7 a_n \cdot T^n / M \quad (\text{A.1})$$

where

c_p is the specific heat capacity, in J · g⁻¹ · K⁻¹;

T is the temperature, in K;

M is the molar mass of corundum $M = 101,961\,3 \text{ g mol}^{-1}$ [3][4];

$a_0, \dots, a_7$ are the coefficients a_0 to a_7 in the temperature intervals according to [Table A.1](#), respectively.

A detailed discussion of the uncertainty of the heat capacity calibration based on α -alumina can be found in Reference [4].

Table A.1 — Polynomial coefficients for calculating the molar heat capacity of corundum (α -Al₂O₃, "Sapphire", molar mass $M = 101,961\,3 \text{ g mol}^{-1}$) [3][4]

	Temperature interval		
	K		
	$70 \leq T < 150$	$150 \leq T < 900$	$900 \leq T < 2\,200$
a_0	$-2,165\,80 \times 10^1$	$-6,877\,17 \times 10^1$	$1,444\,80 \times 10^2$
a_1	$5,093\,80 \times 10^{-1}$	$7,155\,00 \times 10^{-1}$	$-2,385\,35 \times 10^{-2}$
a_2	$1,106\,44 \times 10^4$	$1,456\,85 \times 10^5$	$-8,208\,95 \times 10^6$
a_3	$-1,493\,02 \times 10^{-4}$	$-5,962\,92 \times 10^{-6}$	$2,962\,45 \times 10^{-8}$
a_4	$2,863\,88 \times 10^{-6}$	$1,691\,20 \times 10^{-8}$	$-2,524\,35 \times 10^{-11}$
a_5	$-2,274\,57 \times 10^{-8}$	$-2,210\,44 \times 10^{-11}$	$9,492\,25 \times 10^{-15}$
a_6	$8,538\,76 \times 10^{-11}$	$1,441\,89 \times 10^{-14}$	$-1,702\,62 \times 10^{-18}$
a_7	$-1,256\,07 \times 10^{-13}$	$-3,788\,47 \times 10^{-18}$	$1,158\,82 \times 10^{-22}$