

---

---

**Medical devices — Symbols to be used  
with medical device labels, labelling  
and information to be supplied —**

**Part 1:  
General requirements**

*Dispositifs médicaux — Symboles à utiliser avec les étiquettes,  
l'étiquetage et les informations à fournir relatifs aux dispositifs  
médicaux —*

*Partie 1: Exigences générales*



STANDARDSISO.COM : Click to view the full PDF of ISO 15223-1:2016



**COPYRIGHT PROTECTED DOCUMENT**

© ISO 2016, Published in Switzerland

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office  
Ch. de Blandonnet 8 • CP 401  
CH-1214 Vernier, Geneva, Switzerland  
Tel. +41 22 749 01 11  
Fax +41 22 749 09 47  
copyright@iso.org  
www.iso.org

# Contents

|                                                                                          | Page      |
|------------------------------------------------------------------------------------------|-----------|
| Foreword.....                                                                            | iv        |
| Introduction.....                                                                        | v         |
| <b>1 Scope</b> .....                                                                     | <b>1</b>  |
| <b>2 Normative references</b> .....                                                      | <b>1</b>  |
| <b>3 Terms and definitions</b> .....                                                     | <b>1</b>  |
| <b>4 General requirements</b> .....                                                      | <b>2</b>  |
| 4.1 Proposal of symbols for adoption.....                                                | 2         |
| 4.2 Requirements for usage.....                                                          | 2         |
| 4.3 Other symbols.....                                                                   | 3         |
| <b>5 Symbols</b> .....                                                                   | <b>3</b>  |
| <b>Annex A (informative) Examples</b> .....                                              | <b>20</b> |
| <b>Annex B (informative) Use of general prohibition symbol and negation symbol</b> ..... | <b>24</b> |
| <b>Bibliography</b> .....                                                                | <b>25</b> |

## Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see [www.iso.org/directives](http://www.iso.org/directives)).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see [www.iso.org/patents](http://www.iso.org/patents)).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation on the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT) see the following URL: [www.iso.org/iso/foreword.html](http://www.iso.org/iso/foreword.html).

The committee responsible for this document is ISO/TC 210, *Quality management and corresponding general aspects for medical devices*.

This third edition cancels and replaces the second edition (ISO 15223-1:2012), which has been technically revised with the following principal revisions:

- [Clause 2](#), updated the title of ISO 7000 and added the “date of release” for each of the registered symbols to [Table 1](#);
- symbol 5.1.1, modified the requirement related to the placement of the manufacturer's name and address on IVD labels;
- symbol 5.1.2, modified the requirement related to the placement of name and address of the authorized representative in the European Union on IVD labels;
- symbol 5.4.3, added the information used to indicate an instruction to consult an electronic instructions for use (eIFU);
- symbol 5.4.5, added the reference to ISO 7000, symbol 2725, “Contains or presence of”;
- symbol 5.5.5, modified the description of the symbol and the requirement regarding use with IVD;
- [A.15](#), added the examples of the placement of the eIFU indicator.

A list of all parts in the ISO 15223 series can be found on the ISO website.

NOTE Future symbols intended to appear in this document are to be validated in accordance with ISO 15223-2.

This corrected version of ISO 15223-1:2016 incorporates the following correction:

- in [A.9](#), the graphical symbol of NOTE 2 has been corrected.

## Introduction

This document addresses the presentation of certain items of information that are considered by regulatory authorities to be essential for the safe and proper use of medical devices. As such, the items are required to appear with the medical device in most regulatory domains. The information can be required to appear on the medical device itself, as part of the label, or provided with the medical device.

Many countries require that their own language be used to display textual information with medical devices. At the same time, manufacturers seek to take costs out of labelling by reducing or rationalizing variants. This can cause problems in relation to translation, design and logistics when multiple languages are included on a single label or piece of documentation. For example, users of medical devices labelled in a number of different languages can experience confusion and delay in locating the appropriate language.

This document proposes solutions to these problems through the use of internationally recognized symbols with precisely defined descriptions.

While compiling symbols to be included in this document, ISO/TC 210 recognized the need for systematic methodology for the selection, development and validation of symbols proposed for adoption. This is the subject of ISO 15223-2.

This document is primarily intended to be used by manufacturers of medical devices who market identical products in countries where there are different language requirements for medical device labelling. It can also be of assistance to

- distributors of medical devices or other representatives of manufacturers,
- healthcare providers responsible for training, as well as those being trained,
- those responsible for post-market vigilance,
- healthcare regulatory authorities, testing organizations, certification bodies and other organizations which are responsible for implementing regulations affecting medical devices and which have responsibility for post-market surveillance, and
- consumers or end users of medical devices who draw their supplies from a number of sources and can have varied language capabilities.

STANDARDSISO.COM : Click to view the full PDF of ISO 15223-1:2016

# Medical devices — Symbols to be used with medical device labels, labelling and information to be supplied —

## Part 1: General requirements

### 1 Scope

This document identifies requirements for symbols used in medical device labelling that convey information on the safe and effective use of medical devices. It also lists symbols that satisfy the requirements of this document.

This document is applicable to symbols used in a broad spectrum of medical devices, which are marketed globally and therefore need to meet different regulatory requirements.

These symbols may be used on the medical device itself, on its packaging or in the associated documentation. The requirements of this document are not intended to apply to symbols specified in other standards.

### 2 Normative references

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 7000<sup>1)</sup>, *Graphical symbols for use on equipment — Registered symbols*

ISO 8601, *Data elements and interchange formats — Information interchange — Representation of dates and times*

ISO 15223-2, *Medical devices — Symbols to be used with medical device labels, labelling, and information to be supplied — Part 2: Symbol development, selection and validation*

### 3 Terms and definitions

For the purposes of this document, the terms and definitions given in ISO 14971 and the following apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- IEC Electropedia: available at <http://www.electropedia.org/>
- ISO Online browsing platform: available at <http://www.iso.org/obp>

#### 3.1 characteristic information

information that represents the property or properties of a symbol

1) The collection of ISO 7000 graphical symbols and additional information concerning their use are available at <https://www.iso.org/obp/ui/#search>. Each symbol in the database has a “registration date”. These dates are given in the ISO Registration Number column in Table 1.

### 3.2

#### **description**

normative text which defines the purpose, application and use of the symbol

[SOURCE: IEC 80416-1:2008, 3.2]

### 3.3

#### **label**

written, printed or graphic information provided upon the medical device itself

[SOURCE: GHTE/SG1/N43:2005]

### 3.4

#### **labelling**

information supplied by the manufacturer that is provided for, associated with, or affixed to, a medical device or any of its containers or wrappers

Note 1 to entry: This information relates to the identification, technical *description* (3.2) and use of the medical device, but excludes shipping documents.

Note 2 to entry: Some regional and national regulations refer to “labelling” as “information supplied by the manufacturer”.

### 3.5

#### **symbol used in medical device labelling**

graphical representation appearing on the *label* (3.3) and/or associated documentation of a medical device that communicates *characteristic information* (3.1) without the need for the supplier or receiver of the information to have knowledge of the language of a particular nation or people

Note 1 to entry: The symbol can be an abstract pictorial or a graphical representation, or one that uses familiar objects, including alphanumeric characters.

### 3.6

#### **title**

unique name by which a graphical symbol is identified and spoken of

[SOURCE: IEC 80416-1:2008, 3.9]

## 4 General requirements

### 4.1 Proposal of symbols for adoption

Symbols proposed for adoption in this document shall be validated in accordance with ISO 15223-2.

Any symbol proposed for adoption in this document shall be applicable to a range of medical devices and have global or regional applicability.

### 4.2 Requirements for usage

When risk management shows it to be appropriate for symbols to be used to convey information essential for proper use on the medical device, its packaging or in associated documentation, the symbols given in [Table 1](#) may be used.

Symbols that are registered in ISO 7000 shall comply with the graphical representation in ISO 7000, especially with respect to relative dimensions, including relative line thickness, orientation and the absence or presence of filled or shaded areas.

NOTE 1 ISO and IEC jointly maintain an online database of graphical symbols for use on equipment, which contains the complete set of graphical symbols included in ISO 7000 and IEC 60417 available at <https://www.iso.org/obp/ui/#search>. This online collection shows each graphical symbol and identifies it by a reference number and a title (in English and French). The graphical symbols are available in different formats (e.g. AI, DWG, EPS) and some additional data as applicable is provided. Various search and navigation facilities allow for easy retrieval of graphical symbols.

As part of risk management, the manufacturer should determine the appropriate size for the symbol to be legible for its intended function.

NOTE 2 This document does not specify colours or minimum size for the symbols in [Table 1](#), nor does it specify the relative size of symbols and that of indicated information.

It is important that symbols be used properly. Guidance on appropriate use of the general prohibition symbol and the negation symbol is given in [Annex B](#).

Before symbols are used, the manufacturer shall carry out a risk assessment that indicates that the use of the symbol does not introduce an unacceptable risk.

NOTE 3 Additional information regarding risk assessment can be found in ISO 14971.

Symbols may be used without accompanying text. Where regulations require accompanying text, the title of the symbol given in this document should be considered sufficient. All dates and times presented in association with symbols shall use the conventions set out in ISO 8601.

### 4.3 Other symbols

Other standards specify additional symbols that are applicable to particular kinds or groups of medical devices or to particular situations. Examples of sources for such symbols are identified in the Bibliography. This listing is not exhaustive.

## 5 Symbols

When appropriate, information essential for proper use shall be indicated on the medical device, its packaging, or in the associated documentation by using the corresponding symbols given in [Table 1](#).

A manufacturer may use any appropriate symbol regardless of category.

NOTE [Table 1](#) has been organized into symbol categories for ease of use. The category into which a symbol is grouped does not have any significance as far as usage is concerned. The order of appearance of symbols and the categories in which they are placed are not prioritized. Examples of the use of symbols can be found in [Annex A](#).

Table 1 — Symbols to convey information essential for proper use

| Reference number of symbol | Title of symbol | Description of symbol                                                                                      | Requirements                                                                                                                                                     | Informative notes                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    | Restrictions of use | Additional requirements | ISO 7000 Reg. no. <sup>a</sup> |
|----------------------------|-----------------|------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------|-------------------------|--------------------------------|
| 5.1 Manufacture            |                 |                                                                                                            |                                                                                                                                                                  |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                      |                     |                         |                                |
| 5.1.1                      | Manufacturer    | Indicates the medical device manufacturer, as defined in EU Directives 90/385/EEC, 93/42/EEC and 98/79/EC. | This symbol shall be accompanied by the name and address of the manufacturer (i.e. the person placing the medical device on the market), adjacent to the symbol. | <p>NOTE 1 This symbol is used to indicate information that is required in Europe.<sup>b</sup></p> <p>NOTE 2 The full definition of "manufacturer" is given in EU Directives 90/385/EEC, 93/42/EEC and 98/79/EC.</p> <p>NOTE 3 Guidance on the requirements for EU Directives 90/385/EEC and 93/42/EEC is given in EN 1041.</p> <p>NOTE 4 Guidance on the requirements for EU Directive 98/79/EC is given in ISO 18113-1, ISO 18113-2, ISO 18113-3, ISO 18113-4 and ISO 18113-5.</p> <p>NOTE 5 The date of manufacture, as well as the name and address of the manufacturer, can be combined in one symbol.</p> <p>NOTE 6 The relative size of the symbol and the size of the name and address are not specified.</p> |                     |                         | 3082<br>2011-10-02             |

Table 1 (continued)

| Reference number of symbol                                                                                 | Title of symbol                                     | Description of symbol                                              | Requirements                                                                                                                                 | Informative notes                                                                                                                                                                                                                                                                                                                                                                                                                                        | Restrictions of use | Additional requirements | ISO 7000 Reg. no. <sup>a</sup> |
|------------------------------------------------------------------------------------------------------------|-----------------------------------------------------|--------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------|-------------------------|--------------------------------|
| 5.1.2<br>┌<br><div style="border: 1px solid black; padding: 2px; display: inline-block;">EC REP</div><br>└ | Authorized representative in the European Community | Indicates the authorized representative in the European Community. | This symbol shall be accompanied by the name and address of the authorized representative in the European Community, adjacent to the symbol. | NOTE 1 This symbol is used to indicate information that is required in the European Community.<br>NOTE 2 Guidance on the requirements for EU Directives 90/385/EEC and 93/42/EEC is given in EN 1041.<br>NOTE 3 Guidance on the requirements for EU Directive 98/79/EC is given in ISO 18113-1, ISO 18113-2, ISO 18113-3, ISO 18113-4 and ISO 18113-5.<br>NOTE 4 The relative size of the symbol and the size of the name and address are not specified. |                     |                         | N/A                            |

Table 1 (continued)

| Reference number of symbol                                                                             | Title of symbol     | Description of symbol                                        | Requirements                                                                                                                                                                                                                                                                      | Informative notes                                                                | Restrictions of use | Additional requirements                                                                                                                                                                                                                                                                                                                                                                      | ISO 7000 Reg. no. <sup>a</sup> |
|--------------------------------------------------------------------------------------------------------|---------------------|--------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------|---------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.1.3<br>┌<br><br>└ | Date of manufacture | Indicates the date when the medical device was manufactured. | This symbol shall be accompanied by a date to indicate the date of manufacture. This shall be expressed as in ISO 8601 as four digits for the year and, where appropriate, two digits for the month and two digits for the day. The date shall be located adjacent to the symbol. | NOTE The relative size of the symbol and the size of the date are not specified. |                     | In Europe: <sup>b</sup><br>— the date could be a year, month and day, or year, month and day, as required in the relevant EU Directive;<br>— this symbol may be used to identify the month and year of manufacture for active implantable medical devices, or the year of manufacture for active medical devices where no use-by date is given, as required by the appropriate EU Directive. | 2497<br>2004-01-15             |

Table 1 (continued)

| Reference number of symbol                                                                    | Title of symbol | Description of symbol                                                               | Requirements                                                                                                                                                                                                                                                                                                                                                         | Informative notes                                                                                                                                                                                                                                                                                          | Restrictions of use | Additional requirements                                                                                                                                                                                                                                                                             | ISO 7000 Reg. no. <sup>a</sup> |
|-----------------------------------------------------------------------------------------------|-----------------|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.1.4<br>  | Use-by date     | Indicates the date after which the medical device is not to be used.                | This symbol shall be accompanied by a date to indicate that the medical device should not be used after the end of the year, month or day shown.<br><br>The date shall be expressed as in ISO 8601 as four digits for the year and, where appropriate, two digits for the month and two digits for the day.<br><br>The date shall be located adjacent to the symbol. | NOTE 1 For example, June 2002 is expressed as 2002-06.<br>NOTE 2 The relative size of the symbol and the size of the date are not specified.<br>NOTE 3 Synonym for "use-by date" is "use by".<br>NOTE 4 For some medical devices (e.g. IVDs), this date is only valid when the medical device is unopened. |                     | In Europe: <sup>b</sup><br>— the date could be a year, year and month, or year, month and day, as required by the relevant EU Directive;<br>— this symbol can be used to identify the time limit for implanting an active implantable medical device safely as required by EU Directive 90/385/EEC. | 2607<br>2004-01-15             |
| 5.1.5<br> | Batch code      | Indicates the manufacturer's batch code so that the batch or lot can be identified. | This symbol shall be accompanied by the manufacturer's batch code. The batch code shall be adjacent to the symbol.                                                                                                                                                                                                                                                   | NOTE 1 The relative size of the symbol and the size of the batch code are not specified.<br>NOTE 2 Synonyms for "batch code" are "lot number" and "batch number".                                                                                                                                          |                     |                                                                                                                                                                                                                                                                                                     | 2492<br>2004-01-15             |

Table 1 (continued)

| Reference number of symbol                                                                                  | Title of symbol  | Description of symbol                                                                           | Requirements                                                                                                             | Informative notes                                                                                                                                                                     | Restrictions of use                                                                                                                                                                                                                                                             | Additional requirements | ISO 7000 Reg. no. <sup>a</sup> |
|-------------------------------------------------------------------------------------------------------------|------------------|-------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------|--------------------------------|
| 5.1.6<br>[<br>[<br><br>[ | Catalogue number | Indicates the manufacturer's catalogue number so that the medical device can be identified.     | The manufacturer's catalogue number shall be adjacent to the symbol.                                                     | NOTE 1 The relative size of the symbol and the size of the catalogue number are not specified.<br>NOTE 2 Synonyms for "catalogue number" are "reference number" and "reorder number". | In Europe, <sup>b</sup> the manufacturer's catalogue number shall be placed after or below the symbol and adjacent to it. This symbol may currently be shown without the enclosure; however, it is intended that this option be withdrawn in a future edition of this document. |                         | 2493<br>2004-01-15             |
| 5.1.7<br>[<br>[<br><br>[ | Serial number    | Indicates the manufacturer's serial number so that a specific medical device can be identified. | This symbol shall be accompanied by the manufacturer's serial number. The serial number shall be adjacent to the symbol. | NOTE The relative size of the symbol and the size of the serial number are not specified.                                                                                             | In Europe, <sup>b</sup> the manufacturer's serial number shall be placed after or below the symbol and adjacent to it. This symbol may currently be shown without the enclosure; however, it is intended that this option be withdrawn in a future edition of this document.    |                         | 2498<br>2004-01-15             |

Table 1 (continued)

| Reference number of symbol           | Title of symbol                                | Description of symbol                                                                    | Requirements | Informative notes                                                                                                                                                        | Restrictions of use                                                                                                                                            | Additional requirements | ISO 7000 Reg. no. <sup>a</sup> |
|--------------------------------------|------------------------------------------------|------------------------------------------------------------------------------------------|--------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------|--------------------------------|
| <b>5.2 Sterility</b>                 |                                                |                                                                                          |              |                                                                                                                                                                          |                                                                                                                                                                |                         |                                |
| 5.2.1<br>┌<br><b>STERILE</b><br>└    | Sterile                                        | Indicates a medical device that has been subjected to a sterilization process.           |              | NOTE Use of this symbol precludes the use of symbols 5.2.2 to 5.2.5.                                                                                                     | In Europe, <sup>b</sup> this symbol is restricted to use on terminally sterilized medical devices (EN 556-1:2001, 4.1 applies, including its associated note). |                         | 2499<br>2004-01-15             |
| 5.2.2<br>┌<br><b>STERILE A</b><br>└  | Sterilized using aseptic processing techniques | Indicates a medical device that has been manufactured using accepted aseptic techniques. |              | NOTE 1 Aseptic techniques can include filtration.<br>NOTE 2 Use of this symbol precludes the use of symbol 5.2.1.                                                        |                                                                                                                                                                |                         | 2500<br>2004-01-15             |
| 5.2.3<br>┌<br><b>STERILE EO</b><br>└ | Sterilized using ethylene oxide                | Indicates a medical device that has been sterilized using ethylene oxide.                |              | NOTE Use of this symbol precludes the use of symbol 5.2.1.                                                                                                               | In Europe, <sup>b</sup> this symbol is restricted to use on terminally sterilized medical devices (EN 556-1:2001, 4.1 applies, including its associated note). |                         | 2501<br>2004-01-15             |
| 5.2.4<br>┌<br><b>STERILE R</b><br>└  | Sterilized using irradiation                   | Indicates a medical device that has been sterilized using irradiation.                   |              | NOTE 1 This symbol can be used to indicate that the product has been subjected to irradiation processes.<br>NOTE 2 Use of this symbol precludes the use of symbol 5.2.1. | In Europe, <sup>b</sup> this symbol is restricted to use on terminally sterilized medical devices (EN 556-1:2001, 4.1 applies, including its associated note). |                         | 2502<br>2004-01-15             |

Table 1 (continued)

| Reference number of symbol                                                                                  | Title of symbol                    | Description of symbol                                                                         | Requirements | Informative notes                                                                                                  | Restrictions of use                                                                                                                                            | Additional requirements                                                                                 | ISO 7000 Reg. no. <sup>a</sup> |
|-------------------------------------------------------------------------------------------------------------|------------------------------------|-----------------------------------------------------------------------------------------------|--------------|--------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.2.5<br>「<br>「<br><br>「 | Sterilized using steam or dry heat | Indicates a medical device that has been sterilized using steam or dry heat.                  |              | NOTE Use of this symbol precludes the use of symbol 5.2.1.                                                         | In Europe, <sup>b</sup> this symbol is restricted to use on terminally sterilized medical devices (EN 556-1:2001, 4.1 applies, including its associated note). |                                                                                                         | 2503<br>2004-01-15             |
| 5.2.6<br>「<br><br>「      | Do not resterilize                 | Indicates a medical device that is not to be resterilized.                                    |              |                                                                                                                    |                                                                                                                                                                |                                                                                                         | 2608<br>2004-01-15             |
| 5.2.7<br>「<br><br>「      | Non-sterile                        | Indicates a medical device that has not been subjected to a sterilization process.            |              |                                                                                                                    | This symbol should only be used to distinguish between identical or similar medical devices sold in both sterile and non-sterile conditions.                   |                                                                                                         | 2609<br>2004-01-15             |
| 5.2.8<br>「<br><br>「    | Do not use if package is damaged   | Indicates a medical device that should not be used if the package has been damaged or opened. |              | NOTE This symbol may also mean "Do not use if the product sterile barrier system or its packaging is compromised". |                                                                                                                                                                | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2606<br>2004-01-15             |

Table 1 (continued)

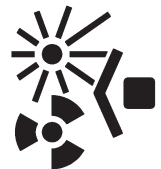
| Reference number of symbol                                                                               | Title of symbol                           | Description of symbol                                                                                                                                                            | Requirements                                                                                                                                                                                         | Informative notes                                                                 | Restrictions of use | Additional requirements                                                                                 | ISO 7000 Reg. no. <sup>a</sup> |
|----------------------------------------------------------------------------------------------------------|-------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|---------------------|---------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.2.9<br>「<br><br>」   | Sterile fluid path                        | Indicates the presence of a sterile fluid path within the medical device in cases when other parts of the medical device, including the exterior, might not be supplied sterile. | The method of sterilization shall be indicated in the empty box, as appropriate. The part of the medical device that is sterile shall be identified in the information supplied by the manufacturer. |                                                                                   |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 3084<br>2011-10-05             |
| <b>5.3 Storage</b>                                                                                       |                                           |                                                                                                                                                                                  |                                                                                                                                                                                                      |                                                                                   |                     |                                                                                                         |                                |
| 5.3.1<br>「<br><br>」   | Fragile, handle with care                 | Indicates a medical device that can be broken or damaged if not handled carefully.                                                                                               |                                                                                                                                                                                                      |                                                                                   |                     |                                                                                                         | 0621<br>2014-06-04             |
| 5.3.2<br>「<br><br>」  | Keep away from sunlight                   | Indicates a medical device that needs protection from light sources.                                                                                                             |                                                                                                                                                                                                      | NOTE This symbol can also mean "Keep away from heat", as referenced in ISO 7000.  |                     |                                                                                                         | 0624<br>2014-06-04             |
| 5.3.3<br>「<br><br>」 | Protect from heat and radioactive sources | Indicates a medical device that needs protection from heat and radioactive sources.                                                                                              |                                                                                                                                                                                                      | NOTE This symbol can also mean "Keep away from sunlight and radioactive sources". |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 0615<br>2004-01-15             |

Table 1 (continued)

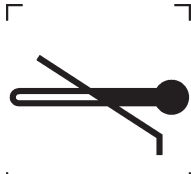
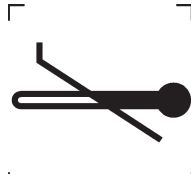
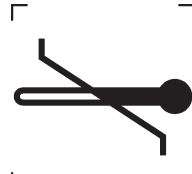
| Reference number of symbol                                                                     | Title of symbol            | Description of symbol                                                                       | Requirements                                                                                                   | Informative notes                                                               | Restrictions of use | Additional requirements | ISO 7000 Reg. no. <sup>a</sup> |
|------------------------------------------------------------------------------------------------|----------------------------|---------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------|---------------------|-------------------------|--------------------------------|
| 5.3.4<br>   | Keep dry                   | Indicates a medical device that needs to be protected from moisture.                        |                                                                                                                | NOTE This symbol can also mean "Keep away from rain" as referenced in ISO 7000. |                     |                         | 0626<br>2014-06-04             |
| 5.3.5<br>   | Lower limit of temperature | Indicates the lower limit of temperature to which the medical device can be safely exposed. | The lower limit of temperature shall be indicated adjacent to the lower horizontal line.                       |                                                                                 |                     |                         | 0534<br>2004-01-15             |
| 5.3.6<br>   | Upper limit of temperature | Indicates the upper limit of temperature to which the medical device can be safely exposed. | The upper limit of temperature shall be indicated adjacent to the upper horizontal line.                       |                                                                                 |                     |                         | 0533<br>2004-01-15             |
| 5.3.7<br> | Temperature limit          | Indicates the temperature limits to which the medical device can be safely exposed.         | The upper and lower limits of temperature shall be indicated adjacent to the upper and lower horizontal lines. |                                                                                 |                     |                         | 0632<br>2014-06-04             |

Table 1 (continued)

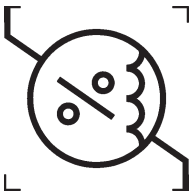
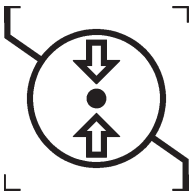
| Reference number of symbol                                                                     | Title of symbol                 | Description of symbol                                                                                              | Requirements                                                                                              | Informative notes                                                       | Restrictions of use | Additional requirements                                                                                 | ISO 7000 Reg. no. <sup>a</sup> |
|------------------------------------------------------------------------------------------------|---------------------------------|--------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------|---------------------|---------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.3.8<br>   | Humidity limitation             | Indicates the range of humidity to which the medical device can be safely exposed.                                 | The humidity limitation shall be indicated adjacent to the upper and lower horizontal lines.              |                                                                         |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2620<br>2004-01-15             |
| 5.3.9<br>   | Atmospheric pressure limitation | Indicates the range of atmospheric pressure to which the medical device can be safely exposed.                     | The atmospheric pressure limitations shall be indicated adjacent to the upper and lower horizontal lines. |                                                                         |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2621<br>2004-01-15             |
| <b>5.4 Safe use</b>                                                                            |                                 |                                                                                                                    |                                                                                                           |                                                                         |                     |                                                                                                         |                                |
| 5.4.1<br>  | Biological risks                | Indicates that there are potential biological risks associated with the medical device.                            | NOTE This symbol is not to be confused with the "Biohazard" sign intended to be used in the workplace.    | See ISO 7010.                                                           |                     |                                                                                                         | 0659<br>2004-01-15             |
| 5.4.2<br> | Do not re-use                   | Indicates a medical device that is intended for one use, or for use on a single patient during a single procedure. |                                                                                                           | NOTE Synonyms for "Do not re-use" are "single use" and "use only once". |                     |                                                                                                         | 1051<br>2004-01-15             |

Table 1 (continued)

| Reference number of symbol                                                                   | Title of symbol              | Description of symbol                                                                                                                                                                                                       | Requirements                                                       | Informative notes                                                                                                                                                                                                                                                                                                          | Restrictions of use                                                                             | Additional requirements                                                                                                                                                                                                                                                                                                                                      | ISO 7000 Reg. no. <sup>a</sup> |
|----------------------------------------------------------------------------------------------|------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.4.3<br> | Consult instructions for use | Indicates the need for the user to consult the instructions for use.                                                                                                                                                        |                                                                    | NOTE 1 Synonym for "Consult instructions for use" is "Consult operating instructions".<br>NOTE 2 Consider the difference between the description of this symbol and that of symbol 5.4.4.                                                                                                                                  |                                                                                                 | When used to indicate an instruction to consult an electronic instructions for use (eIFU), this symbol is accompanied by an eIFU indicator. This indicator may represent the manufacturer's eIFU website or any other appropriate indication on the use of eIFU. The indicator may be placed either alongside, beneath or surrounding the symbol (see A.15). | 1641<br>2004-01-15             |
| 5.4.4<br> | Caution                      | Indicates the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself. | The symbol variant ISO 7000: symbol 0434B ("Caution") may be used. | NOTE 1 Consider the difference between the description of this symbol and that of symbol 5.4.3.<br>NOTE 2 This symbol is essentially a cautionary symbol and is used to highlight the fact that there are specific warnings or precautions associated with the medical device, which are not otherwise found on the label. | This symbol is not to be confused with the "Caution" sign intended to be used in the workplace. |                                                                                                                                                                                                                                                                                                                                                              | 0434A<br>2004-01-15            |

Table 1 (continued)

| Reference number of symbol                                                                               | Title of symbol                              | Description of symbol                                                                                                                                               | Requirements | Informative notes                                                                                               | Restrictions of use                                                                                                                                              | Additional requirements                                                                                                           | ISO 7000 Reg. no. <sup>a</sup>                                                                                                     |
|----------------------------------------------------------------------------------------------------------|----------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------|-----------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------|
| 5.4.5<br>「<br><br>」   | Contains or presence of natural rubber latex | Indicates the presence of natural rubber or dry natural rubber latex as a material of construction with in the medical device or the packaging of a medical device. |              | NOTE This symbol is intended to warn those people who may have allergic reactions to certain proteins in latex. | This symbol should not be used for medical devices containing “synthetic rubber”.                                                                                | In Europe, b this symbol shall be explained in the information supplied by the manufacturer.                                      | Application of ISO 7000, symbol 2725<br>2005-09-08                                                                                 |
| 5.5 IVD-specific                                                                                         |                                              |                                                                                                                                                                     |              |                                                                                                                 |                                                                                                                                                                  |                                                                                                                                   |                                                                                                                                    |
| 5.5.1<br>「<br><br>」   | In vitro diagnostic medical device           | Indicates a medical device that is intended to be used as an <i>in vitro</i> diagnostic medical device.                                                             |              |                                                                                                                 | This symbol should only be used to identify <i>in vitro</i> diagnostic medical devices and not to specify that the medical device is for “ <i>in vitro</i> use”. | In Europe, b this symbol is only used to identify <i>in vitro</i> diagnostic medical devices as defined in EU Directive 98/79/EC. | N/A                                                                                                                                |
| 5.5.2<br>「<br><br>」  | Control                                      | Indicates a control material that is intended to verify the performance characteristics of another medical device.                                                  |              | NOTE For negative controls, use symbol 5.5.3 and for positive controls, use symbol 5.5.4.                       |                                                                                                                                                                  |                                                                                                                                   | N/A<br>NOTE ISO 7000-2494 was withdrawn by ISO/TC 145/SC 3. Letters and words are not registered as graphical symbols in ISO 7000. |
| 5.5.3<br>「<br><br>」 | Negative control                             | Indicates a control material that is intended to verify the results in the expected negative range.                                                                 |              |                                                                                                                 |                                                                                                                                                                  |                                                                                                                                   | 2495<br>2004-01-15                                                                                                                 |

Table 1 (continued)

| Reference number of symbol                                                                                       | Title of symbol                   | Description of symbol                                                                               | Requirements                                                                                | Informative notes                                                                | Restrictions of use | Additional requirements | ISO 7000 Reg. no. <sup>a</sup>                     |
|------------------------------------------------------------------------------------------------------------------|-----------------------------------|-----------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------|---------------------|-------------------------|----------------------------------------------------|
| 5.5.4<br>┌<br>└<br><br>┌<br>└ | Positive control                  | Indicates a control material that is intended to verify the results in the expected positive range. |                                                                                             |                                                                                  |                     |                         | 2496<br>2004-01-15                                 |
| 5.5.5<br>┌<br>└<br><br>┌<br>└ | Contains sufficient for <n> tests | Indicates the total number of IVD tests that can be performed with the IVD.                         | The number of tests that can be performed with the IVD shall appear adjacent to the symbol. | NOTE The relative size of the symbol and the number of tests performed can vary. |                     |                         | Application of ISO 7000, symbol 0518<br>2004-01-15 |

Table 1 (continued)

| Reference number of symbol                                                                               | Title of symbol                     | Description of symbol                                                                                                                                                                          | Requirements | Informative notes                                                                                                                                                                                                                                                | Restrictions of use                                                                                                                                                                                                                                                                                                                                                                            | Additional requirements                                                                                 | ISO 7000 Reg. no. <sup>a</sup>                     |
|----------------------------------------------------------------------------------------------------------|-------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------|----------------------------------------------------|
| 5.5.6<br>[<br><br>]   | For IVD performance evaluation only | Indicates an IVD device that is intended to be used only for evaluating its performance characteristics before it is placed on the market for medical diagnostic use.                          |              | NOTE 1 A synonym is "IVD for investigational use only".<br>NOTE 2 A medical device that is for IVD performance evaluation only is not intended to be used for an <i>in vitro</i> diagnostic examination for medical purposes (i.e. to yield diagnostic results). | This symbol shall not appear jointly on the label or in the labelling of an IVD device bearing the symbol<br><br>[<br><br>]<br><br>which means that the medical device is an <i>in vitro</i> diagnostic medical device intended by the manufacturer to be used for an <i>in vitro</i> diagnostic examination. |                                                                                                         | Application of ISO 7000, symbol 3083<br>2011-10-03 |
| 5.6 Transfusion/infusion                                                                                 |                                     |                                                                                                                                                                                                |              |                                                                                                                                                                                                                                                                  |                                                                                                                                                                                                                                                                                                                                                                                                |                                                                                                         |                                                    |
| 5.6.1<br>[<br><br>] | Sampling site                       | Indicates a medical device or blood processing application that includes a system dedicated to the collection of samples of a given substance stored in the medical device or blood container. |              | NOTE This is not to be associated with a site on a patient where samples are taken.                                                                                                                                                                              |                                                                                                                                                                                                                                                                                                                                                                                                | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2715<br>2005-09-08                                 |

Table 1 (continued)

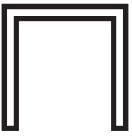
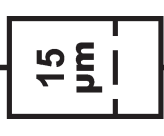
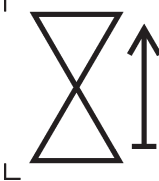
| Reference number of symbol                                                                     | Title of symbol              | Description of symbol                                                                                                       | Requirements | Informative notes                                                                                                                                       | Restrictions of use | Additional requirements                                                                                 | ISO 7000 Reg. no. <sup>a</sup> |
|------------------------------------------------------------------------------------------------|------------------------------|-----------------------------------------------------------------------------------------------------------------------------|--------------|---------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------|---------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.6.2<br>   | Fluid path                   | Indicates the presence of a fluid path.                                                                                     |              |                                                                                                                                                         |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2722<br>2005-09-08             |
| 5.6.3<br>   | Non-pyrogenic                | Indicates a medical device that is non-pyrogenic.                                                                           |              |                                                                                                                                                         |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2724<br>2005-09-08             |
| 5.6.4<br>   | Drops per millilitre         | Indicates the number of drops per millilitre.                                                                               |              | NOTE The number of drops per millilitre is specified; 20 is shown as an example and will be replaced by the appropriate number of drops per millilitre. |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2726<br>2005-09-08             |
| 5.6.5<br> | Liquid filter with pore size | Indicates an infusion or transfusion system of the medical device that contains a filter of a particular nominal pore size. |              | NOTE The nominal pore size of the filter is specified; 15 is shown as an example and will be replaced by the appropriate pore size.                     |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2727<br>2005-09-08             |

Table 1 (continued)

| Reference number of symbol                                                                                                                                                                                                                                                                                                                                           | Title of symbol | Description of symbol                                                           | Requirements | Informative notes                                                                                                 | Restrictions of use | Additional requirements                                                                                 | ISO 7000 Reg. no. <sup>a</sup> |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------|---------------------------------------------------------------------------------|--------------|-------------------------------------------------------------------------------------------------------------------|---------------------|---------------------------------------------------------------------------------------------------------|--------------------------------|
| 5.6.6<br>                                                                                                                                                                                                                                                                         | One-way valve   | Indicates a medical device with a valve that allows flow in only one direction. |              | NOTE It is important for the user to know that the flow is only possible in one direction and cannot be reversed. |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2728<br>2005-09-08             |
| 5.7 Others                                                                                                                                                                                                                                                                                                                                                           |                 |                                                                                 |              |                                                                                                                   |                     |                                                                                                         |                                |
| 5.7.1<br>                                                                                                                                                                                                                                                                         | Patient number  | Indicates a unique number associated with an individual patient.                |              |                                                                                                                   |                     | In Europe, <sup>b</sup> this symbol shall be explained in the information supplied by the manufacturer. | 2610<br>2004-01-15             |
| <sup>a</sup> Each symbol in the ISO/IEC symbols database (available at <a href="https://www.iso.org/obp/ui">https://www.iso.org/obp/ui</a> ) has a "registration date". This date is shown below the Registration Number.<br><sup>b</sup> At the moment, only countries applying the principles laid down in the EU Directives have this requirement or restriction. |                 |                                                                                 |              |                                                                                                                   |                     |                                                                                                         |                                |