
**Fine ceramics (advanced ceramics,
advanced technical ceramics) —
Determination of bacterial reduction
rate by semiconducting photocatalytic
materials under indoor lighting
environment — Semi-dry method
for estimating antibacterial activity
on the actual environmental bacteria
contamination surface**

STANDARDSISO.COM : Click to view the full PDF of ISO 22551:2020



STANDARDSISO.COM : Click to view the full PDF of ISO 22551:2020



COPYRIGHT PROTECTED DOCUMENT

© ISO 2020

All rights reserved. Unless otherwise specified, or required in the context of its implementation, no part of this publication may be reproduced or utilized otherwise in any form or by any means, electronic or mechanical, including photocopying, or posting on the internet or an intranet, without prior written permission. Permission can be requested from either ISO at the address below or ISO's member body in the country of the requester.

ISO copyright office
CP 401 • Ch. de Blandonnet 8
CH-1214 Vernier, Geneva
Phone: +41 22 749 01 11
Fax: +41 22 749 09 47
Email: copyright@iso.org
Website: www.iso.org

Published in Switzerland

Contents

	Page
Foreword	iv
Introduction	v
1 Scope	1
2 Normative references	1
3 Terms and definitions	1
4 Principle	2
5 Materials	2
5.1 Bacteria strains and preparation for tests	2
5.1.1 Bacteria strains	2
5.1.2 Subculture of bacteria strains	3
5.2 Chemicals and implements	3
5.2.1 Dispersion medium for <i>Staphylococcus epidermidis</i>	3
5.2.2 Dispersion medium for <i>Escherichia coli</i>	3
5.2.3 1/5-concentration nutrient broth medium (1/5 NB)	3
5.2.4 1/2-concentration luris-bertani broth medium (1/2 LB)	3
5.2.5 Nutrient agar medium (NA medium)	4
5.2.6 Soybean-casein digest broth with lecithin and polysorbate 80 (SCDLP medium)	4
6 Apparatus	4
7 Test piece	5
8 Procedure	5
8.1 Preparation of test bacterial solution	5
8.1.1 Culturing of test bacteria	5
8.1.2 Preparation of test bacterial solution	5
8.2 Application of test bacterial solution	5
8.2.1 Measurement of test piece mass	5
8.2.2 Application of test bacterial solution	5
8.2.3 Direct application method	5
8.3 Indoor lighting condition	6
8.3.1 Measurement of luminance and preparation of setting position of test pieces	6
8.3.2 Light exposure time	6
8.3.3 Test pieces inoculated with test bacterial solution under light condition	6
8.3.4 Test pieces inoculated with test bacterial solution under dark condition	6
8.4 Measurement of number of living bacteria	7
8.4.1 Washout of inoculated bacteria	7
8.4.2 Measurement of bacterial concentration in washout solution	7
8.4.3 Measurement of bacterial concentration in test bacterial solution	7
9 Calculation	8
9.1 General	8
9.2 Test requirement fulfilment validation	8
9.3 Calculation of reduction rate of living bacteria	9
10 Test report	10
Bibliography	12

Foreword

ISO (the International Organization for Standardization) is a worldwide federation of national standards bodies (ISO member bodies). The work of preparing International Standards is normally carried out through ISO technical committees. Each member body interested in a subject for which a technical committee has been established has the right to be represented on that committee. International organizations, governmental and non-governmental, in liaison with ISO, also take part in the work. ISO collaborates closely with the International Electrotechnical Commission (IEC) on all matters of electrotechnical standardization.

The procedures used to develop this document and those intended for its further maintenance are described in the ISO/IEC Directives, Part 1. In particular, the different approval criteria needed for the different types of ISO documents should be noted. This document was drafted in accordance with the editorial rules of the ISO/IEC Directives, Part 2 (see www.iso.org/directives).

Attention is drawn to the possibility that some of the elements of this document may be the subject of patent rights. ISO shall not be held responsible for identifying any or all such patent rights. Details of any patent rights identified during the development of the document will be in the Introduction and/or on the ISO list of patent declarations received (see www.iso.org/patents).

Any trade name used in this document is information given for the convenience of users and does not constitute an endorsement.

For an explanation of the voluntary nature of standards, the meaning of ISO specific terms and expressions related to conformity assessment, as well as information about ISO's adherence to the World Trade Organization (WTO) principles in the Technical Barriers to Trade (TBT), see www.iso.org/iso/foreword.html.

This document was prepared by Technical Committee ISO/TC 206, *Fine ceramics*.

Any feedback or questions on this document should be directed to the user's national standards body. A complete listing of these bodies can be found at www.iso.org/members.html.

Introduction

This document applies to testing the antibacterial activity of indoor-light-active photocatalytic ceramics and other materials, produced by either coating or mixing of a light-active photocatalyst in a simulated indoor lighting environment. The International Standard for testing the antibacterial activity of photocatalytic materials has been published as ISO 27447 and the International Standard for testing the antibacterial activity of indoor-light-active photocatalytic materials has been published as ISO 17094.

Recently, it was discovered that antibacterial activities by indoor-light-active photocatalyst showed different results between the ISO 17094 testing in a laboratory and the evaluation in a real environment. Thus, the ISO 17094 testing method should be complemented by a testing method that simulates an indoor lighting environment.

ISO 17094 is one of the best methods for measuring the potential of antibacterial activity using photocatalysts. This document can measure the reduction rate of living bacteria by testing with impure substances similar to an indoor lighting environment.

[STANDARDSISO.COM](https://standardsiso.com) : Click to view the full PDF of ISO 22551:2020

Fine ceramics (advanced ceramics, advanced technical ceramics) — Determination of bacterial reduction rate by semiconducting photocatalytic materials under indoor lighting environment — Semi-dry method for estimating antibacterial activity on the actual environmental bacteria contamination surface

WARNING — Handling and manipulation of microorganisms that are potentially hazardous requires a high degree of technical competence. Only personnel trained in microbiological techniques should carry out tests.

1 Scope

This document establishes a test method for determining the antibacterial activity of materials containing an indoor-light-active photocatalytic material on the surface. The antibacterial reduction rate is determined by measuring the survival of bacteria after illumination with indoor light. This test assumes a surface with high potential of possible person contact with bacteria. This test is designed to evaluate the suppression of contact infection of bacteria using an indoor-light-active photocatalytic material under indoor lighting environment.

It is intended for use with different kinds of indoor-light-active photocatalytic materials used in construction, for example, flat sheets, board or plate shapes, which are the basic forms of materials for various applications. It is not applicable to powder, granular, or porous indoor-light-active photocatalytic materials, as well as cloths or textiles.

It is applicable to indoor-light-active photocatalytic materials produced for antibacterial application. Other types of indoor-light-active photocatalytic materials applications, i.e. decomposition of water contaminants, self-cleaning, antifogging, and air purification, are non-applicable by this method.

2 Normative references

The following documents are referred to in the text in such a way that some or all of their content constitutes requirements of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 14605, *Fine ceramics (advanced ceramics, advanced technical ceramics) — Light source for testing semiconducting photocatalytic materials used under indoor lighting environment*

ISO 17094:2014, *Fine ceramics (advanced ceramics, advanced technical ceramics) — Test method for antibacterial activity of semiconducting photocatalytic materials under indoor lighting environment*

3 Terms and definitions

For the purposes of this document, the following terms and definitions apply.

ISO and IEC maintain terminological databases for use in standardization at the following addresses:

- ISO Online browsing platform: available at <https://www.iso.org/obp>
- IEC Electropedia: available at <http://www.electropedia.org/>

3.1

photocatalyst

substance that performs one or more functions based on oxidization and reduction reactions under photoirradiation, including decomposition and removal of air and water contaminants, deodorization, antibacterial, antifungal, self-cleaning and antifogging actions

3.2

indoor-light-active photocatalyst

photocatalyst (3.1) that functions under illumination with artificial light used for general lighting purposes

3.3

indoor lighting environment

illumination with artificial light source(s) used for general lighting purposes and excluding sunlight

3.4

indoor-light-active photocatalytic material

material in which or on which the *indoor-light-active photocatalyst* (3.2) is added by, for example, coating, impregnation or mixing

3.5

reduction rate of living bacteria

reduction of the viable bacteria counts on the surface of *indoor-light-active photocatalytic materials* (3.4) versus that of non-treated materials after treatment under illumination or in a dark condition for a designated time, expressed as a proportion of the bacteria in percentage per unit of inoculating mass

3.6

dispersion medium

medium used for the preparation of a test bacterial solution, comprising a chemical component, which simulates human sebum

4 Principle

The testing method is used to obtain the antibacterial activity of indoor-light-active photocatalytic materials under indoor lighting condition by contact of a test piece with bacteria suspended in a dispersion medium simulating human sebum. This method is applicable to materials with smooth surfaces, such as flat sheet, board or plate-shaped materials.

The surface of the test piece is uniformly inoculated with bacterial suspension. The inoculated test piece is placed on a Petri dish and a glass plate is placed on top of the Petri dish. The Petri dish containing the test piece is exposed to light. After exposure, the test bacteria are washed out of the test piece. This washout suspension is measured by the viable bacterial count method.

5 Materials

5.1 Bacteria strains and preparation for tests

5.1.1 Bacteria strains

The bacterial strains to be used in the test shall be the same as or equivalent to those described in [Table 1](#) and supplied by an entry registered under the World Federation for Culture Collections or the Japan Society for Microbial Resources and Systematics.

Table 1 — Bacteria strains to be used in the test

Bacteria species	WDCM code
<i>Staphylococcus epidermidis</i>	WDCM 00132
<i>Escherichia coli</i>	WDCM 00012

5.1.2 Subculture of bacteria strains

Testing is performed in a biosafety level testing room equipped for the required type of microorganisms to be used. A nutrient agar slant medium (NA slant medium) is inoculated with *Staphylococcus epidermidis* or *Escherichia coli* and streak cultured for 24 h to 48 h at $(37 \pm 1)^\circ\text{C}$. The cultured bacteria can be stored at 5°C to 10°C for one month. Using this material as a stock strain, the cultured bacteria can be subcultured up to 10 times under the same conditions.

NOTE 1 For long-term storage such as cryopreservation, appropriate measures are taken.

NOTE 2 If the activity of used bacteria is maintained, agar plates can be used.

5.2 Chemicals and implements

The chemicals, materials and implements used in this testing are as follows, unless otherwise designated. The commercial media of the same components described below can be used.

5.2.1 Dispersion medium for *Staphylococcus epidermidis*

Add 49 ml of sterile water and 50 g of weighed triol-type polypropylene glycol (alpha, alpha', alpha"-1,2,3-propanetriyltris [omega-hydroxy-poly[oxy(methyl-1,2-ethanediyl)]], average molecular weight 300 (PPG 300)) to a flask with a magnetic stirring bar. Immediately after addition, stir with the bar. Then, add 1 g of polyethylene glycol [average molecular weight 500 000 (PEG 500K)] under continuous stirring, and stir for 1 h. Sterilize in an autoclave at $(121 \pm 2)^\circ\text{C}$ for at least 15 min. After dissolving, use sodium hydroxide or hydrochloric acid to adjust the solution pH to $(7,1 \pm 0,2)$ at 25°C . Adjust to 100 ml using sterile water. Sterilize in an autoclave at $(121 \pm 2)^\circ\text{C}$ for at least 15 min and store at room temperature. Do not use dispersion medium made more than one month previously.

5.2.2 Dispersion medium for *Escherichia coli*

Add 30 ml of sterile water and 0,1 g of weighed mucin to a flask with a magnetic stirring bar. Immediately after addition, stir with the bar. Then, add 0,2 g of sodium alginate under continuous stirring. After dissolving, add 20 ml of glycerol. Store at room temperature. Do not use dispersion medium made more than one month previously.

5.2.3 1/5-concentration nutrient broth medium (1/5 NB)

For 100 ml of sterile water, take 0,3 g of meat extract, 1,0 g of peptone and 0,5 g of sodium chloride, put them into a flask and dissolve them thoroughly. When the contents are thoroughly dissolved, use a solution of sodium hydroxide or hydrochloric acid to bring the pH to $(7,1 \pm 0,1)$ at 25°C . Dilute this solution fivefold using sterile water, and adjust the pH to $(7,0 \pm 0,2)$ at 25°C using sodium hydroxide or hydrochloric acid. Sterilize in an autoclave at $(121 \pm 2)^\circ\text{C}$ for at least 15 min. After preparation, if 1/5 NB is not used immediately, store it at 5°C to 10°C . Do not use 1/5 NB made more than one month previously.

5.2.4 1/2-concentraion luris-bertani broth medium (1/2 LB)

For 100 ml of sterile water, take 0,5 g yeast extract, 1 g tryptone and 1 g sodium chloride, put them into a flask and dissolve them thoroughly. When the contents are thoroughly dissolved, use a solution of sodium hydroxide or hydrochloric acid to bring the pH to $(7,1 \pm 0,1)$ at 25°C . Dilute this solution twofold

using sterile water, and adjust the pH to $(7,0 \pm 0,2)$ at 25 °C using sodium hydroxide or hydrochloric acid. Sterilize in an autoclave at (121 ± 2) °C for at least 15 min. After preparation, if 1/2 LB is not used immediately, store it at 5 °C to 10 °C. Do not use 1/2 LB made more than one month previously.

5.2.5 Nutrient agar medium (NA medium)


The preparation of the NA medium shall follow ISO 17094:2014, 6.2.3.

5.2.6 Soybean-casein digest broth with lecithin and polysorbate 80 (SCDLP medium)

The preparation of the SCDLP medium shall follow ISO 17094:2014, 6.2.4.

6 Apparatus

The testing equipment enables examination of antibacterial activity using an indoor-light-active photocatalytic material by providing indoor light illumination to activate the indoor-light-active photocatalyst. It consists of a light source and a chamber containing the test piece. A schematic representation of the testing system is shown in [Figure 1](#).

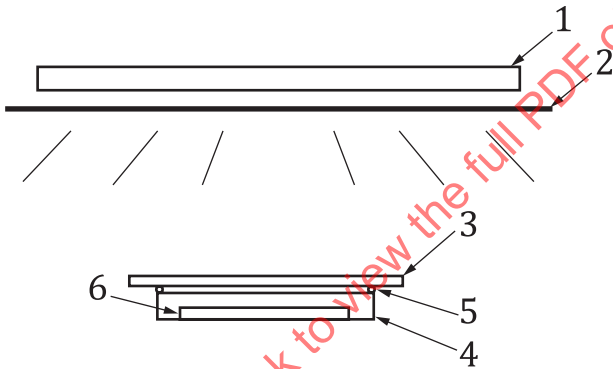

Key

- | | | | |
|---|-------------------------|---|------------|
| 1 | light source | 4 | Petri dish |
| 2 | UV sharp cut-off filter | 5 | spacer |
| 3 | glass plate | 6 | test piece |

Figure 1 — Schematic diagram of the test equipment

6.1 Light source, as specified in ISO 14605, for indoor lighting condition. A halophosphate fluorescent lamp with a correlated colour temperature of 3 800 K to 4 500 K shall be used for the light source. When a halophosphate fluorescent lamp is not available, a three-band fluorescent lamp with a correlated colour temperature of 3 800 K to 4 500 K and a colour rendering index (Ra) higher than 80 can be used as a substitute.

6.2 UV sharp cut-off filter, as specified in ISO 14605, used under UV cut-off condition (condition B). The UV sharp cut-off filter shall be mounted immediately below the lamp. Condition B is under 380 nm cut-off condition. Type B UV sharp cut-off-filter is specified in ISO 14605.

6.3 Moisture preservation glass plate, as specified in ISO 17094:2014, 7.3, to prevent contamination by floating bacteria.

6.4 Illuminance meter, as specified in ISO 14605.

7 Test piece

Cut a (50 ± 2) mm $\times$ (50 ± 2) mm flat portion of the material. The thickness of the material should not exceed 10 mm. Use this as the standardized shaped test piece. A film-like material may be attached to non-treated glass to maintain its strength. Prepare nine pieces of non-treated test pieces and six pieces of indoor-light-active photocatalytic treated test pieces. When non-treated test pieces cannot be provided, a material as similar as possible to that of the indoor-light-active photocatalytic treated test pieces can be used. Take great care to avoid microbial contamination and cross-contamination among test pieces. If the surface of the test piece is contaminated by organic material, the test piece may be UV irradiated for up to 24 h using approximately 1,0 mW/cm², or wiped clean with ethanol. If a non-treated test piece cannot be provided, it is necessary to use a material that would produce a similar bacteria recovery rate to that of an indoor-light-active photocatalytic treated test piece.

8 Procedure

8.1 Preparation of test bacterial solution

8.1.1 Culturing of test bacteria

Transfer bacteria (5.1.2) to the NA medium (slant or plate) using a platinum loop and incubate at (37 ± 1) °C for 24 h to 48 h. Transfer the incubate bacteria to a new NA medium (slant or plate) and incubate at (37 ± 1) °C for 24 h to 48 h. This cultivated bacteria shall be used as test bacteria.

8.1.2 Preparation of test bacterial solution

Scrape the test bacteria cultured in 8.1.1 from the NA medium using a platinum loop, suspend the material in a glass reagent tube containing 5 ml of 1/5 NB for *Staphylococcus epidermidis*, 10 ml of 1/2 LB for *Escherichia coli*, and agitate thoroughly. The obtained solution shall be used as the original bacterial solution. Add the dispersion medium described in 5.2.1 (for *Staphylococcus epidermidis*) or 5.2.2 (for *Escherichia coli*) to the glass reagent tube. Then add 1/100 (for *Staphylococcus epidermidis*) or 1/10 (for *Escherichia coli*) of the amount of the original bacterial solution and agitate thoroughly. Adjust the test bacterial solution to a concentration of $5,0 \times 10^2$ cfu/mg to $2,0 \times 10^3$ cfu/mg. After adjustment, use immediately. If the test bacterial solution is not to be used immediately, store it on ice and use within 2 h. In all cases, use glass reagent tubes for preparation of the test bacterial solution. *Staphylococcus epidermidis* adheres readily to the surface of plastic materials.

NOTE The high viscosity of the dispersion medium impedes blending of the dispersion medium and the original bacterial solution. One way to achieve a nearly homogenous test bacterial solution is to repeat the process three times or more (inverting the tube five times and mixing the tube by vortex mixer for 5 s).

8.2 Application of test bacterial solution

8.2.1 Measurement of test piece mass

Measure the mass of all test pieces to a precision of 0,1 mg and record.


8.2.2 Application of test bacterial solution

Apply $(2,5 \pm 0,5)$ mg of the test bacterial solution as prepared in 8.1.2 to test pieces (nine non-treated test pieces and six indoor-light-active photocatalytic treated test pieces) uniformly. 8.2.3 provides an example of the application method. However, other similar application methods may be employed.

8.2.3 Direct application method

Transfer 0,2 ml of the bacterial solution to the tip of a sterilized cellulose non-woven folded cloth as shown in Figure 2. Prepare a piece of aluminium foil by wiping clean with alcohol and inoculate an area of approximately 15 cm $\times$ 15 cm to eliminate excess bacterial solution. Then apply the solution on the

test piece surface completely and uniformly. After the application is complete, wipe off any excess test bacterial solution uniformly with a newly folded cellulose non-woven cloth. Calculate the mass of the applied solution by measuring the mass of the test piece on which the solution has been applied. Repeat either application or removal by wiping until the application mass is $(2,5 \pm 0,5)$ mg. Place the inoculated test piece into a Petri dish and cover with a glass plate.

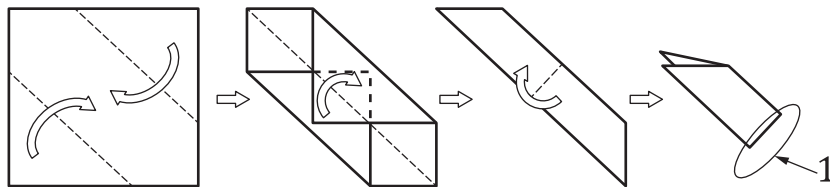

Key

1 section to be impregnated with test bacterial solution

Figure 2 — Folding of cellulose non-woven cloth

Three test pieces comprise one group. Limit the time of the bacterial solution application for one group to 5 min or less. Then proceed to expose to light or keep in a dark place, immediately.

8.3 Indoor lighting condition

8.3.1 Measurement of luminance and preparation of setting position of test pieces

Set the illuminance meter on the base of the illumination apparatus and place a glass plate used for testing on top of the sensor. Adjust the intensity of the lamp to give an illuminance of $1\,000\text{ lx} \pm 50\text{ lx}$ at the surface of the test pieces. Before measuring luminance, the light source shall have been lighting up preliminarily for 15 min or longer to stabilize the luminance. The testing ambient temperature and relative humidity shall be set within the range of $(25 \pm 3)^\circ\text{C}$ and 50 % to 70 %, respectively, to simulate the environment where the product may be used.

NOTE This illuminance range can be altered between $100\text{ lx} \pm 5\text{ lx}$ and $3\,000\text{ lx} \pm 150\text{ lx}$ to take into account the real conditions where the indoor-light-active photocatalyst product is applied.

8.3.2 Light exposure time

Expose to light the Petri dishes containing the test pieces (three non-treated and three indoor-light-active photocatalytic treated test pieces) with bacterial inoculation for 4 h.

NOTE This exposure time can be shortened to as little as 2 h to take into account the real conditions where the indoor-light-active photocatalyst product is applied.

8.3.3 Test pieces inoculated with test bacterial solution under light condition

Place the test pieces (three non-treated and three indoor-light-active photocatalytic treated test pieces) with inoculated test bacterial solution on Petri dishes and place a glass plate on top of each Petri dish. Place Petri dishes containing the test pieces at the position determined in 8.3.1 and expose to light for the duration designated in 8.3.2.

8.3.4 Test pieces inoculated with test bacterial solution under dark condition

Place test pieces (three non-treated and three indoor-light-active photocatalytic treated test pieces) with inoculated test bacterial solution on Petri dishes and place a glass plate on top of each Petri dish. Keep the Petri dishes containing the test pieces in a dark place for the same time as instructed in 8.3.3.

8.4 Measurement of number of living bacteria

8.4.1 Washout of inoculated bacteria

Put each tested piece (non-treated and indoor-light-active photocatalytic treated test pieces) on a new Petri dish using sterilized tweezers. Add 10 ml of SCDLP medium (see 5.2.6) kept at 37 °C using a sterilized pipette. Place the dish on a horizontal shaker and shake for 5 min at a rate of 90 cycles per min. Another method may be used if a similar yield is assured.

Use a Petri dish large enough that the surface of a test piece is completely submerged when 10 ml of SCDLP medium is added.

8.4.2 Measurement of bacterial concentration in washout solution

Take 1,0 ml of washout solution (see 8.4.1) with a sterilized pipette, add to a reagent tube containing $(9 \pm 0,1)$ ml of SCDLP medium and mix thoroughly. Take three Petri dishes and place in each 1,0 ml of each solution (washout solution or diluted solution). Add 15 ml to 20 ml of NA medium (see 5.2.5) kept at 45 °C to 48 °C to each Petri dish. After mixing well, cover and let it stand for 15 min at room temperature. After the medium solidifies, place the Petri dishes upside down and culture in an incubator at (37 ± 1) °C for 36 h to 48 h. Count colony numbers in the series of Petri dishes with 30 colonies to 300 colonies. The bacterial concentration in washout solution is calculated using Formula (1) to two significant digits.

$$P = Z \times D_F \quad (1)$$

where

- P is the bacterial concentration (cfu/ml);
- Z is the average number of colonies in the three Petri dishes (cfu);
- D_F is the dilution factor.

If the number of viable bacteria is less than 30 in the Petri dishes with 1,0 ml of washout solution, count all and the average is used for the calculation. If the number of viable bacteria is less than one in the Petri dishes with 1,0 ml of the washout solution, the number of viable bacteria is designated as "<1" and Z is replaced with "1" for the calculation.

8.4.3 Measurement of bacterial concentration in test bacterial solution

Measure the mass of a (glass or plastic) reagent tube containing $(9 \pm 0,1)$ ml of SCDLP medium, then add approximately 1,0 g of the test bacterial solution to the reagent tube. Measure the mass again and record the mass of the added test bacterial solution [designate mass as M_A (g)], then mix thoroughly. Take 1 ml and add to a reagent tube containing $(9 \pm 0,1)$ ml of SCDLP medium. Repeat this process to prepare a tenfold serial dilution. Take two Petri dishes and place in each one 1,0 ml of each solution (diluted solution). Add 15 ml to 20 ml of NA medium (see 5.2.5) kept at 45 °C to 48 °C to each Petri dish. After mixing well, cover and let it stand for 15 min at room temperature. After the medium solidifies, place the Petri dishes upside down and culture in an incubator at (37 ± 1) °C for 36 h to 48 h. Count colony numbers in the series of Petri dishes with 30 colonies to 300 colonies. The bacterial concentration in test bacterial solution is calculated using Formula (2) to two significant digits.

$$C = Z_2 \times D_F / 1\,000 / M_A \quad (2)$$

where

C is the bacterial concentration in test bacterial solution (cfu/mg);

Z_2 is the average number of colonies on the two Petri dishes (cfu);

D_F is the dilution factor;

M_A is the mass of added test bacterial solution (g).

NOTE The specific gravity of the dispersion medium and water is the same.

9 Calculation

9.1 General

The test results are calculated as follows. The calculated values should be rounded to the second decimal place in accordance with ISO 80000-1.

9.2 Test requirement fulfilment validation

A test is considered valid if it fulfils all of the following four conditions from a) to d). If one or more of these conditions are not fulfilled, the test is considered non-valid and shall be performed again.

The number of viable bacteria per unit inoculating mass on each of the test pieces (N) shall be calculated using [Formula \(3\)](#).

$$N = P \times V / M_1 \quad (3)$$

where

N is the number of viable bacteria per unit inoculating mass on each test piece (cfu/mg);

P is the concentration of viable bacteria (cfu/ml) obtained in [8.4.2](#);

V is the volume of SCDLP medium used for washout (ml);

M_1 is the mass of inoculate on each test piece (mg).

- a) The bacterial concentration in the test bacterial solution (C) shall be within the range described in [Formula \(4\)](#).

$$5,0 \times 10^2 \leq C \leq 2,0 \times 10^3 \quad (4)$$

where C is the bacterial concentration in the test bacterial solution (cfu/mg).

- b) Distribution of the average viable bacteria count per unit inoculating mass of each set of the three non-treated test pieces under illumination with light or kept in a dark condition after inoculation shall be within the range described in [Formula \(5\)](#).

$$(N_{\max} - N_{\min}) / N_{\text{mn}} \leq 0,6 \quad (5)$$

where

- $N_{\max}$ is the maximum value of viable bacteria per unit inoculating mass on the three individual test pieces;
- $N_{\min}$ is the minimum value of viable bacteria per unit inoculating mass on the three individual test pieces;
- N_{mn} is the average value of viable bacteria per unit inoculating mass on the three individual test pieces.

- c) The average value of viable bacteria per unit inoculating mass among the three non-treated test pieces immediately after inoculation shall be within the range described in [Formula \(6\)](#).

$$0,5 \times C \leq A \leq 1,3 \times C \quad (6)$$

where

- C is the bacterial concentration in the test bacterial solution (cfu/mg);
- A is the average viable bacteria count per unit inoculating mass among the three non-treated test pieces immediately after inoculation (cfu/mg).
- d) The ratio of the average viable bacteria count per unit inoculating mass among the three non-treated test pieces immediately after inoculation (A) to the average viable bacteria count per unit inoculating mass among the three non-treated test pieces exposed to light illumination or kept in a dark condition (B_L or B_D) shall be within the range described in [Formulae \(7\)](#) and [\(8\)](#), respectively.

$$0,5 \leq (B_L / A) \leq 1,3 \quad (7)$$

$$0,5 \leq (B_D / A) \leq 1,3 \quad (8)$$

where

- B_L is the average viable bacteria count per unit inoculating mass among the three non-treated test pieces exposed to light illumination (cfu/mg);
- B_D is the average viable bacteria count per unit inoculating mass among the three non-treated test pieces kept in a dark condition (cfu/mg);
- A is the average viable bacteria count per unit inoculating mass among the three non-treated test pieces immediately after inoculation (cfu/mg).

9.3 Calculation of reduction rate of living bacteria

The reduction rate of living bacteria on indoor-light-active photocatalytic treated test pieces exposed to light illumination (R_L) and the reduction rate of living bacteria on indoor-light-active photocatalytic treated test pieces kept in a dark condition (R_D) shall be calculated with [Formulae \(9\)](#) and [\(10\)](#), respectively.

$$R_L = 100 \times (B_L - C_L) / B_L \quad (9)$$

where

R_L is the reduction rate of living bacteria on the indoor-light-active photocatalytic treated test pieces exposed to light illumination (%);

B_L is the average viable bacteria count per unit inoculating mass among the three non-treated test pieces exposed to light illumination (cfu/mg);

C_L is the average viable bacteria count per unit inoculating mass among the three indoor-light-active photocatalytic treated test pieces exposed to light illumination (cfu/mg).

$$R_D = 100 \times (B_D - C_D) / B_D \quad (10)$$

where

R_D is the reduction rate of living bacteria on the indoor-light-active photocatalytic treated test pieces kept in a dark condition (%);

B_D is the average viable bacteria count per unit inoculating mass among the three non-treated test pieces kept in a dark condition (cfu/mg);

C_D is the average viable bacteria count per unit inoculating mass among the three indoor-light-active photocatalytic treated test pieces kept in a dark condition (cfu/mg).

10 Test report

The test report shall include the following information:

- a) a reference to this document, i.e. ISO 22551:2020;
- b) description of the type, size, shape and thickness of the indoor-light-active photocatalytic treated test pieces and non-treated test pieces;
- c) description of conditions of pre-exposure when applied;
- d) type of test bacteria, bacteria strain number;
- e) manufacturer of fluorescent lamp, product number;
- f) type of UV sharp cut-off filter (manufacturer, product number);
- g) manufacturer of illuminance meter, product number;
- h) lighting conditions including luminance and light exposure duration;
- i) type and size of glass plate;
- j) composition of dispersion medium;
- k) bacterial concentration in test bacterial solution;
- l) inoculating method;
- m) washout method;
- n) values obtained by testing (R_L , R_D , A , B_L , C_L , B_D , C_D);