



SURFACE VEHICLE STANDARD

J2020™

OCT2022

Issued 1989-06
Revised 2022-10

Superseding J2020 APR2016

Accelerated Exposure of Automotive Exterior Materials Using a Fluorescent UV and Condensation Apparatus

RATIONALE

Updated references including CIE 241; clarifications on black panel thermometers; specimen repositioning guidelines; and various grammatical and wording changes.

1. SCOPE

- 1.1 This test method specifies the operating conditions for a fluorescent ultraviolet (UV) and condensation apparatus used for the accelerated exposure of various automotive exterior components.
- 1.2 Specimen preparation, test duration, and performance evaluation procedures are addressed by each automotive manufacturer's material specifications.
- 1.3 This SAE Standard may involve hazardous materials, operations, and equipment. This document does not purport to address all of the safety problems associated with its use.

It is the responsibility of whoever uses this document to consult and establish appropriate and health practices and determine the applicability of regulatory limitations prior to use.

1.4 Significance and Use

This test method is designed to simulate extreme environmental conditions encountered on the outside of an automobile due to sunlight, heat, and moisture, and to provide an acceleration of exposure for the purpose of predicting the performance of exterior automotive materials.

- 1.5 Unless otherwise specified, all dimensions reported in this method are nominal.

SAE Executive Standards Committee Rules provide that: "This report is published by SAE to advance the state of technical and engineering sciences. The use of this report is entirely voluntary, and its applicability and suitability for any particular use, including any patent infringement arising therefrom, is the sole responsibility of the user."

SAE reviews each technical report at least every five years at which time it may be revised, reaffirmed, stabilized, or cancelled. SAE invites your written comments and suggestions.

Copyright © 2022 SAE International

All rights reserved. No part of this publication may be reproduced, stored in a retrieval system or transmitted, in any form or by any means, electronic, mechanical, photocopying, recording, or otherwise, without the prior written permission of SAE.

TO PLACE A DOCUMENT ORDER: Tel: 877-606-7323 (inside USA and Canada)
Tel: +1 724-776-4970 (outside USA)
Fax: 724-776-0790
Email: CustomerService@sae.org
http://www.sae.org

SAE WEB ADDRESS:

For more information on this standard, visit
https://www.sae.org/standards/content/J2020_202210/

2. REFERENCES

2.1 Applicable Documents

The following publications form a part of this specification to the extent specified herein. Unless otherwise indicated, the latest issue of SAE publications shall apply.

2.1.1 ASTM Publications

Available from ASTM International, 100 Barr Harbor Drive, P.O. Box C700, West Conshohocken, PA 19428-2959, Tel: 610-832-9585, www.astm.org.

- | | |
|-----------|--|
| ASTM E207 | Standard Test Method for Thermal EMF Test of Single Thermoelement Materials by Comparison with a Reference Thermoelement of Similar EMF-Temperature Properties |
| ASTM E220 | Method for Calibration of Thermocouples by Comparison Techniques |
| ASTM G113 | Standard Terminology Relating to Natural and Artificial Weathering Tests of Nonmetallic Materials |
| ASTM G154 | Standard Practice for Operating Fluorescent Light Apparatus for UV Exposure of Nonmetallic Materials |

2.1.2 CIE Publications

Available from CIE Central Bureau, Babenbergerstrasse 9/9A, 1010 Vienna, Austria, Tel: +43 1 714 31 87, www.cie.co.at.

CIE Pub No. 241 Recommended Reference Solar Spectra for Industrial Applications

2.1.3 ISO Publications

Copies of these documents are available online at <http://webstore.ansi.org/>.

- | | |
|----------|--|
| ISO 9370 | Plastics - Instrumental Determination of Radiant Exposure in Weathering Tests - General Guidance and Basic Test Method |
|----------|--|

2.2 Related Publications

The following publications are provided for information purposes only and are not a required part of this SAE Technical Report.

2.2.1 ASTM Publications

Available from ASTM International, 100 Barr Harbor Drive, P.O. Box C700, West Conshohocken, PA 19428-2959, Tel: 610-832-9585, www.astm.org.

- | | |
|-----------|--|
| ASTM G151 | Standard Practice for Exposing Nonmetallic Materials in Accelerated Test Devices that Use Laboratory Light Sources |
|-----------|--|

2.2.2 ISO Publications

Copies of these documents are available online at <http://webstore.ansi.org/>.

- | | |
|-------------|--|
| ISO 4892-1 | Plastics - Methods of Exposure to Laboratory Light Sources - Part 1: General Guidance |
| ISO 16474-1 | Plastics - Instrumental Determination of Radiant Exposure in Weathering Tests - General Guidance and Basic Test Method |
| ISO 16474-3 | Paints and Varnishes - Methods of Exposure to Laboratory Light Sources - Part 3: Fluorescent UV Lamps |

3. DEFINITIONS

The terminology found in ASTM G113 may be appropriate to this procedure.

3.1 BLACK PANEL THERMOMETER, n

A temperature measuring device used to control an artificial weathering device and to provide an estimate of the maximum temperature of specimens exposed to a radiant energy source.

3.2 IRRADIANCE, n

The radiant power per unit area incident on a receiver, typically reported in watts per square meter (W/m^2).

3.3 IRRADIANCE, SPECTRAL, n

The distribution of irradiance as a function of wavelength ($\text{W}/\text{m}^2/\text{nm}$).

3.4 SPECTRAL POWER DISTRIBUTION (SPD), n

The absolute or relative radiant power emitted by a source, or incident upon a receiver as a function of wavelength.

3.5 FLUORESCENT UV LAMP, n

A lamp in which the irradiance from a low-pressure mercury arc is transformed to a higher wavelength UV by a phosphor. The spectral power distribution of a fluorescent lamp is determined by the emission spectrum of the phosphor and the UV transmittance of the glass tube.

3.6 APPARATUS

A more complete description of the apparatus may be found in ASTM G154.

3.6.1 Test chamber, constructed of corrosion-resistant materials enclosing eight fluorescent UV lamps, a heated water pan, test specimen racks, and provisions for controlling and indicating operating times and temperatures.

3.7 Lamps shall be rapid-start, medium bi-pin fluorescent UV type with a length of 1200 mm, and a nominal rating of 40 W.

3.7.1 The lamps shall be one of the following types:

a. UVB lamps with a peak emission of 313 nm as described in Table 3 of ASTM G154. A representative spectral power distribution of this lamp is shown in Figure 1.

NOTE: In devices that do not automatically control irradiance, "FS-40" or "F40UVB" type lamps have historically been used. In devices that automatically control irradiance, UVB-313 lamps may be used.

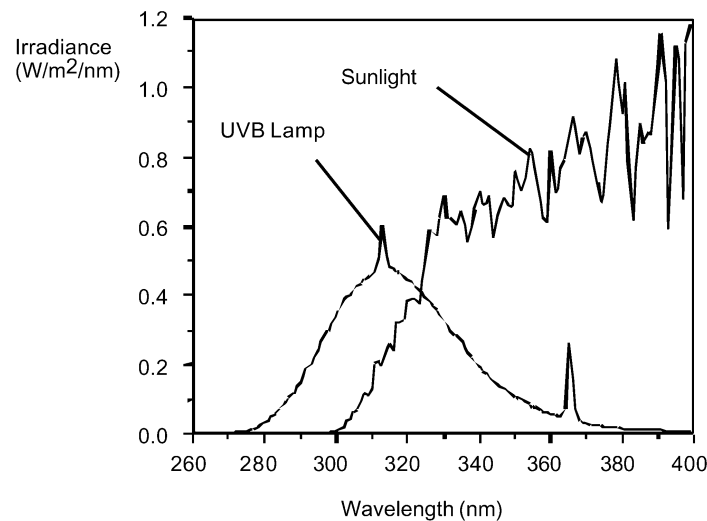


Figure 1 - UVB lamps representative spectral power distribution

3.7.2 Other fluorescent UV lamps meeting the size and electrical characteristics in 3.7 may be used by prior agreement, provided that the lamp and spectral power distribution are reported in conformance with the report section. Use of lamps other than those specified in 3.7 may result in significant differences in test results.

- a. UVA lamps with a peak emission of 340 nm as described in Table 1 of ASTM G154. A representative spectral power distribution of this lamp is shown in Figure 2.

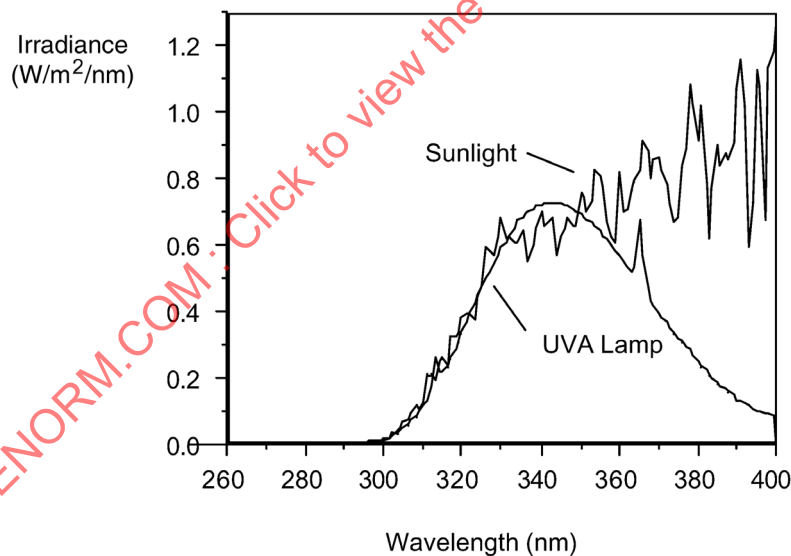


Figure 2 - UVA lamps representative spectral power distribution

3.8 LAMP SPACING AND ARRANGEMENT

The lamps shall be mounted in two banks of four lamps each. The lamps in each bank shall be mounted in a flat plane on 70 mm centers.

- 3.8.1 Lamp manufacturing techniques typically result in a slight gradient in light output from one end of a lamp to the other end. To compensate for this, see Figure 3.

3.9 TEST SPECIMENS

The test specimens shall be mounted in stationary racks with the plane of the surface parallel to the plane of the lamps at a distance of 50 mm from the nearest surface of the lamps.

3.10 CONDENSATION MECHANISM

Water vapor shall be generated by heating a water pan extending under the entire specimen area and containing a minimum water depth of 20 mm. Specimen racks and the test specimens themselves shall constitute the side walls of the chamber. The back side of the specimens shall be exposed to cooling effects of ambient room air. The resulting heat transfer causes water to condense on the test surface.

3.10.1 WATER SUPPLY

Supply water with an automatic control to regulate the level in the water pan shall be provided. Distilled, deionized, or potable tap water are equally acceptable for purposes of the test, since the condensation process itself distills water onto the test surface.

NOTE: Use of distilled or deionized water is recommended. Using tap water may necessitate more frequent cleaning of water pan.

3.11 CYCLE TIME

A continuously operating cycle timer shall be provided to program the selected cycle of UV periods and condensation periods.

3.12 SPECIMEN TEMPERATURE MEASUREMENT AND CONTROL

Specimen temperature shall be measured and controlled by a black panel thermometer. This consists of a black-coated metal panel that absorbs all wavelength uniformly, with a thermally sensitive element firmly attached to the center of the exposed surface. Typical dimensions are 75 × 100 × 2.5 mm thick. The temperature sensor shall be accurate to ± 2 °C through the temperature range defined in Section 5. The temperature indicator shall be located outside the test chamber. The black panel thermometer shall be positioned so that the sensor is subject to the same conditions as the specimens.

During UV exposure, the selected equilibrium temperature shall be maintained within ± 3 °C of set-point temperature by supplying heated air to the test chamber. During condensation exposure, the selected equilibrium temperature shall be maintained within ± 3 °C of set-point temperature by heating the water in the water pan. The UV and condensation temperature controls shall be independent of each other.

3.13 TEST CHAMBER LOCATION

3.13.1 The apparatus shall be located in an area maintained between 18 °C and 27 °C. Control of ambient temperature is particularly critical for proper operation of apparatus stacked one above the other, because the heat generated from one unit could interfere with the operation of adjacent units.

3.13.2 The apparatus shall be located at least 300 mm from walls or other apparatus. Nearby heat sources, such as ovens or heated test apparatus, shall be avoided or shielded.

3.13.3 The room where the apparatus is located shall be ventilated to remove the heat and moisture produced and to maintain the temperatures specified in 5.5.

3.14 To ensure repeatability of test results, maintain and calibrate the apparatus to the manufacturer's instructions.

4. TEST SPECIMENS

- 4.1 Specimen size shall be either 50 × 75 mm, 75 × 100 mm, 75 × 150 mm, 100 × 150 mm, or as agreed upon by concerned parties.

NOTE: Because of the geometry of some specimen holders, the previous sizes may cause some problems if, after inspection, they are replaced in a different orientation than the original. For example, improper replacement of 75 × 150 mm panel into a commonly used holder can cause the end 1/3 to be shielded (i.e., exposed less than the center 1/3).

- 4.1.1 For simultaneous exposures of original and repair coatings on one test specimen, 100 × 150 mm specimens have been found useful. In this situation, each coating should cover one half of the specimen and be 50 × 150 mm area.
- 4.2 Replicate specimens are desirable to provide a record of degradation at different time intervals. Retention of unexposed specimens in a controlled environment is recommended as it is difficult to mask a specimen to prevent exposure to condensation.
- 4.3 For specimens of insulating materials, such as plastics or foams, maximum specimen thickness should be 30 mm to allow adequate heat transfer for condensation. Report any thickness greater than this in the test report.
- 4.4 To provide rigidity, flexible specimens should be attached to a backing panel. This panel can, for example, be made of aluminum, 0.635 mm thick, 3003 H14 alloy.
- 4.5 Holes in specimens larger than 2 mm and any openings larger than 1 mm around irregularly shaped specimens shall be sealed to prevent loss of water vapor. Porous specimens, such as textiles, shall be backed with a vapor barrier such as metal.

5. PROCEDURE

- 5.1 Mount the test specimens in the specimen racks with the test surfaces facing the lamp. When the test specimens do not completely fill the racks, the empty spaces must be filled with panels to maintain the test conditions within the chamber. Gaps between specimen holders greater than 1 mm can affect test conditions and are not permitted.
- 5.2 The test specimens shall be exposed within an area 210 mm in height by 900 mm wide on each side of the apparatus located as shown in Figure 3.

NOTE: Specimens in the extreme left and right side of the tester (positions 1 and 13 in Figure 4) experience a lower irradiance than other specimens. While these positions do meet the irradiance requirements in ASTM G151 when repositioning is used, it is recommended that these positions are excluded when test and control specimens do not completely fill the specimen racks.

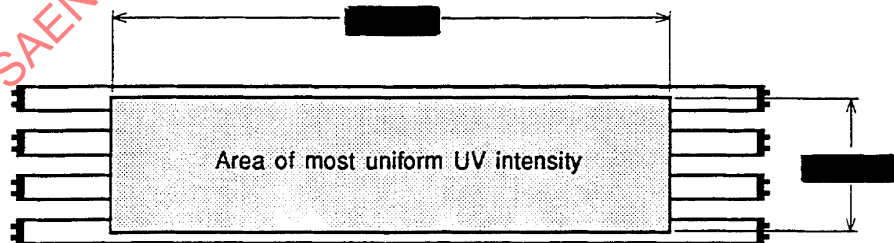


Figure 3 - Area of most uniform UV intensity

- 5.3 Set condensation cooling timer to 15 minutes, if applicable.
- 5.4 Initiate all exposures at the beginning of the condensation cycle.

- 5.5 Program the cycle timer to achieve the following test conditions: 8 hours UV light exposure at 70 °C, alternating with 4 hours condensation exposure at 50 °C.
- 5.5.1 Operating continuously, repeating the cycle, except for servicing the instrument and inspection of specimens. Inspect specimens weekly during the condensation cycle for evidence of condensation.
- 5.6 In order to minimize any effect from temperature or UV light variation, specimens shall be repositioned periodically. As a minimum, specimens must be rotated horizontally by (a) moving the two extreme right-hand specimen holders to the far left of the exposure area, and (b) sliding the remaining specimen holders to the right (see Figure 4). Repositioning must be performed a minimum of four times per exposure test. For exposures with a duration longer than approximately 1 month, it is recommended to reposition at least weekly.

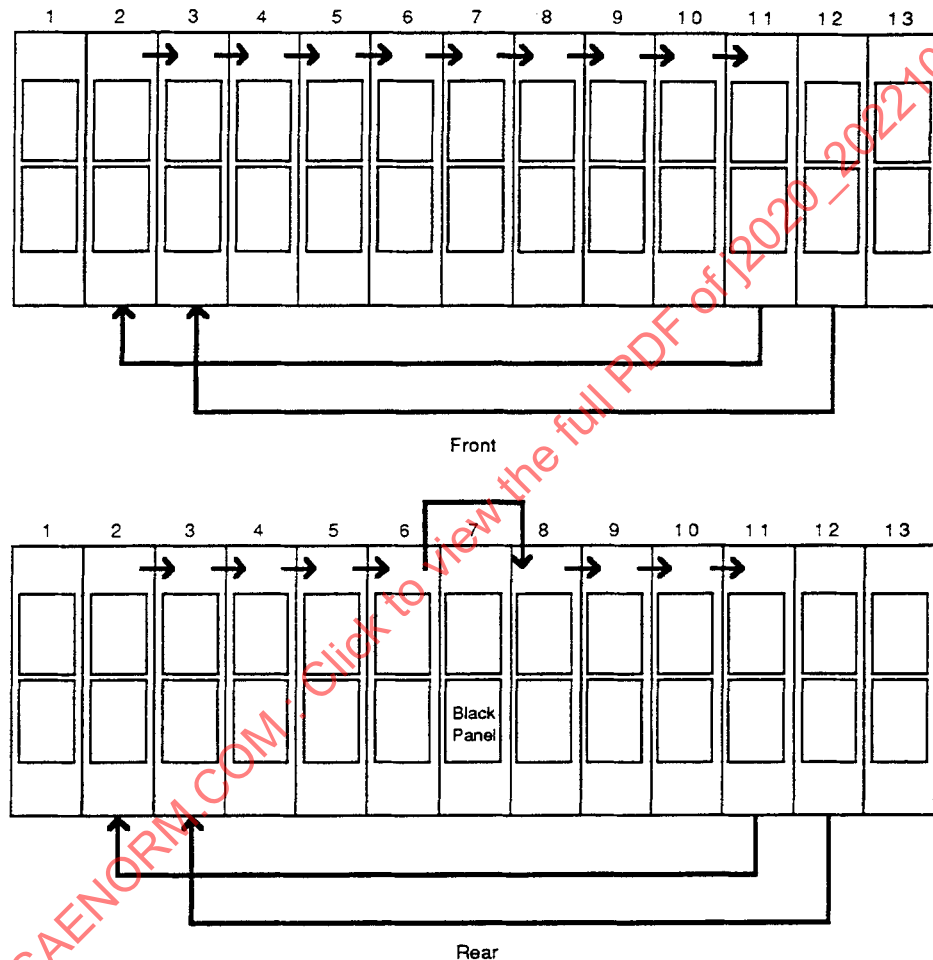


Figure 4 - Horizontal specimen holder rotation
(see 5.2 for guidance on use of end positions)

5.7 UV Irradiance

- 5.7.1 UV irradiance shall be maintained within a specified range by one of the methods shown in Appendix A.
- 5.7.2 For UVB lamps in chambers with irradiance control, the irradiance level shall be 0.48 W/m²/nm at 310 nm, with allowable operational fluctuation of ±0.02 W/m²/nm. For chambers without irradiance control, irradiance shall be periodically measured and maintained between 0.41 W/m²/nm and 0.55 W/m²/nm at 310 nm by lamp replacement and/or rotation (see A.2.3). If the irradiance deviates beyond this tolerance, terminate the test until the cause of the deviation has been determined and corrected.
- 5.7.2.1 Other irradiance levels can be used by mutual agreement between contracting parties.

5.7.3 For UVA lamps, the irradiance level shall be mutually agreed upon by contracting parties.

NOTE: For devices that automatically control irradiance, target irradiance levels between the range of 0.67 to 0.89 W/m²/nm at 340 nm are commonly used. CIE Publication No. 241, Table A2, states that noon summer sunlight, air mass 1.0, is 0.76 W/m²/nm at 340 nm.

5.8 Expose specimens for the required time. Refer to material/manufacturer specification for specific requirements.

6. REPORT

6.1 The report may be prepared using Figure 5, if appropriate. The report shall include the following:

6.1.1 Laboratory.

6.1.2 Material.

6.1.3 Test method.

6.1.4 The supplier and model of fluorescent UV/Condensation apparatus.

6.1.5 The name of the supplier of the fluorescent UV lamp. The manufacturer's designation for the lamp, the lot number or date code, and the wavelength (nm) at which peak emission occurs. For example: FS-40, C6 lot, 313 nm.

6.1.6 Cycle of UV exposure time and temperature, condensation time and temperature. For example: 8 hours UV/70 °C, 4 hours condensation/50 °C.

6.1.7 Total hours exposure time.

6.1.8 Results of UV irradiance measured and the intervals at which it was measured during the exposure test as specified in Appendix A.

6.1.9 Any special conditions of test.

SAENORM.COM : Click to view the full PDF of J2020_202210

**ACCELERATED EXPOSURE OF AUTOMOTIVE EXTERIOR MATERIALS
USING A FLUORESCENT UV AND CONDENSATION APPARATUS**

Material Test Date: _____

Material Identification: _____

Total Hours Exposure: _____

Test Cycle_____ hours UV at _____ °C
_____ hours Condensation at _____ °C**UV Light Source**

Lamp Type: _____ Supplier: _____

Lot Number: _____ Peak Emission: _____

Total Time Meter at Start of Test: _____ Total Time Meter at End of Test: _____

Frequency of Sample Rotation: _____

Apparatus Model: _____

Special Test Conditions: _____
_____Remarks: _____

Name/Signature: _____

Date: _____

Attach measurements of irradiance. Include all individual measurements and the average irradiance throughout the exposure test.*Figure 5 - Exposure report form***7. NOTES****7.1 Revision Indicator**

A change bar (|) located in the left margin is for the convenience of the user in locating areas where technical revisions, not editorial changes, have been made to the previous issue of this document. An (R) symbol to the left of the document title indicates a complete revision of the document, including technical revisions. Change bars and (R) are not used in original publications, nor in documents that contain editorial changes only.

APPENDIX A - PERIODIC MAINTENANCE - NON-MANDATORY

A.1 MAINTENANCE AND CALIBRATION OF TEST APPARATUS

- A.1.1 The test apparatus requires periodic maintenance to maintain uniform exposure conditions. Perform required maintenance and calibration in accordance with manufacturer's instructions.

A.2 MAINTENANCE OF IRRADIANCE

- A.2.1 Maintenance of proper UV irradiance can be achieved by any of the two methods specified as follows. Regardless of the method used for regular monitoring and control, the system should be measured every 6 months using the procedure specified in Appendix B.

A.2.2 Method 1, Automatic Feedback-Loop Method

The system consists of a programmable controller that continuously monitors the UV intensity. A feedback loop system maintains the programmed irradiance level by adjusting power to the UV lamps.

- A.2.2.1 Every 500 hours $\pm$ 25 hours of lamp operation, measure light output and calibrate the integrated (on-board) irradiance measurement and control system with a calibration radiometer as recommended by the manufacturer.

- A.2.2.2 Replace all of the lamps in the unit when they can no longer maintain the irradiance setpoint.

- A.2.2.3 Using the procedure in B.1.3, re-calibrate the calibration radiometer at a minimum every 12 months.

A.2.3 Method 2, Manual Irradiance Measurement Method

This system uses constant output ballasts which do not compensate for variability in light output caused by lamp aging, temperature, or lot-to-lot variability.

NOTE: Many fluorescent lamps age significantly with extended use. It is recommended to follow the apparatus manufacturer's instructions on the procedure necessary to maintain desired irradiance.

- A.2.3.1 At a minimum, operators should measure lamp irradiance every 2 weeks (336 total hours) to assure that the irradiance is within the specified range. If the light intensity is below $0.48 \pm 15\%$ W/m²/nm at 310 nm (0.41 to 0.55 W/m²/nm at 310 nm), replace one lamp in each bank and rotate the others as shown in Figure A1.

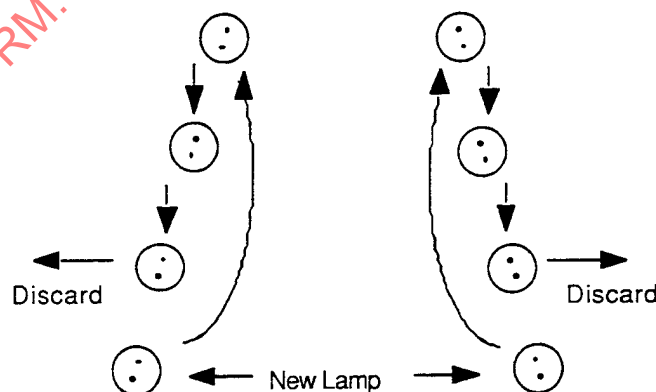


Figure A1 - Lamp rotation/replacement procedure

- A.2.3.2 Using the procedure in B.1.3, re-calibrate the calibration radiometer at a minimum every 12 months.